

# Tutorial: FIB preparation of lamellae on chips for in-situ TEM and other techniques

J. Reuteler<sup>1\*</sup>

<sup>1</sup> ETH Zurich, ScopeM (Scientific Center for Optical and Electron Microscopy), 8093 Zürich, Switzerland

\*email: [joakim.reuteler@scopem.ethz.ch](mailto:joakim.reuteler@scopem.ethz.ch)

This tutorial targets a wide range of FIB operators, from beginners to long time expert users across fields, interested in preparation of lamellae on chips (LoC). Starting with a recap of the in-situ lift-out technique for preparing a TEM lamella mounted to a half-grid (Omnigrid), the terminology used throughout the tutorial is defined. Important concepts for achieving high quality TEM lamellae for different purposes are summarized.

There is an increasing demand for thin samples mounted to a specific chip, typically MEMS. Thus the core part of this tutorial discusses several protocols for preparation of lamellae on chips. Examples will cover the preparation of samples for in-situ TEM observations, synchrotron-based x-ray techniques under local magnetic field as well as for heating experiments, and a lamella precisely mounted on an optical waveguide. Furthermore, different routes for dealing with unconventional geometries are discussed.

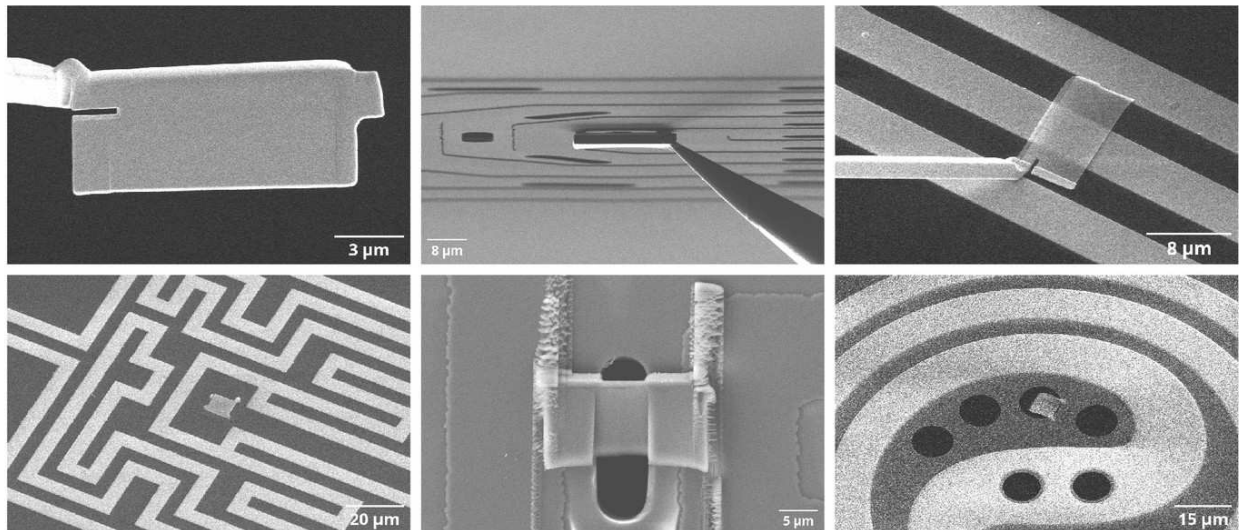


Fig. 1: Snapshots taken while preparing various lamellae on chip.

# Tungsten based SQUID Nanofabrication by means of Focused Ion Beam Induced Deposition

F. Sigloch<sup>1</sup>, S. Soraya<sup>1,2</sup>, P. Orús<sup>1,2</sup> and J. M. De. Teresa<sup>1,2</sup>

<sup>1</sup> Instituto de Nanociencia y Materiales de Aragón (INMA), Universidad de Zaragoza-CSIC, 50009 Zaragoza, Spain

<sup>2</sup> Laboratorio de Microscopía Avanzadas (LMA), Universidad de Zaragoza, 50018 Zaragoza, Spain

\* corresponding author email: [fsigloch@unizar.es](mailto:fsigloch@unizar.es)

Direct write techniques such as Focused Electron/Ion Beam Induced Deposition (FEBID/FIBID) constitute versatile, resist-free techniques for the fabrication of nanostructures, offering an alternative to conventional methods such as Optical or Electron Beam Lithography.

The deposition of the commercially available precursor gas  $W(CO)_6$  with  $Ga^+$  ions results in a film of W-C with well-established superconducting properties [1]. Planar deposits exhibit a critical temperature of  $T_c = 4\text{ K} - 5\text{ K}$ , an upper critical magnetic field of  $B_{c,2} = 7\text{ T} - 8.5\text{ T}$  and a critical current density of  $J_c = 0.001\text{ MA/cm}^2 - 0.01\text{ MA/cm}^2$  [2]. The London penetration depth is reported to be  $\lambda_L = 850\text{ nm}$  and the superconducting coherence length as  $\xi = 6\text{ nm} - 9\text{ nm}$  [2]. With  $Ga^+$  FIBID and the  $W(CO)_6$  precursor, superconducting nanostructures with a linewidth of 50 nm are feasible to fabricate with high precision and reproducibility.

In this work we present the fabrication of nanoscale Superconducting Quantum Interference Devices (nanoSQUIDs) of W-C by means of  $Ga^+$  FIBID. The SQUID loop is formed by a 50 nm thick and 200 nm wide film with a rectangular inner loop of  $300 \times 700\text{ nm}^2$ . The Josephson Junctions (JJs) are formed by 300 nm long constrictions with a cross-sectional area of down to  $50 \times 50\text{ nm}^2$ . The SQUIDs obtained show a critical temperature of up to  $T_c = 4.3\text{ K}$  and a critical current of up to  $I_c = 8.5\text{ }\mu\text{A}$ . The normal state resistance of the JJs is  $R_N = 496\text{ }\Omega$ . Upon variation of the external magnetic field we observe periodic oscillations in the critical current,  $I_c$ , and upon injection of a constant bias current  $I_b \sim I_c$  the voltage dropping across the structure displays a sinusoidal dependence on the external magnetic field. The transfer coefficient is remarkably high, with up to  $V_\phi = 1301\text{ }\mu\text{A}/\phi_0$  due to the high normal-state resistance of the JJs [3]. Recent efforts towards improving the nanoSQUIDs properties will be introduced.

[1] P. Orús, F. Sigloch, S. Sangiao and J. M. De Teresa. *Superconducting Materials and Devices Grown by Focused Ion and Electron Beam Induced Deposition*, Nanomaterials 12, (2022) 1367

[2] P. Orús, R. Córdoba, J. M. De Teresa. *Nanofabrication – Nanolithography techniques and their application*, IOP Publishing (2020), ch. 5.

[3] F. Sigloch, S. Sangiao, P. Orús, J. M. De Teresa. *Large output voltage to magnetic flux change in nanoSQUIDs based on direct-write Focused Ion Beam Induced Deposition technique*, arXiv [Preprint] (2022), [arXiv:2203.05278].

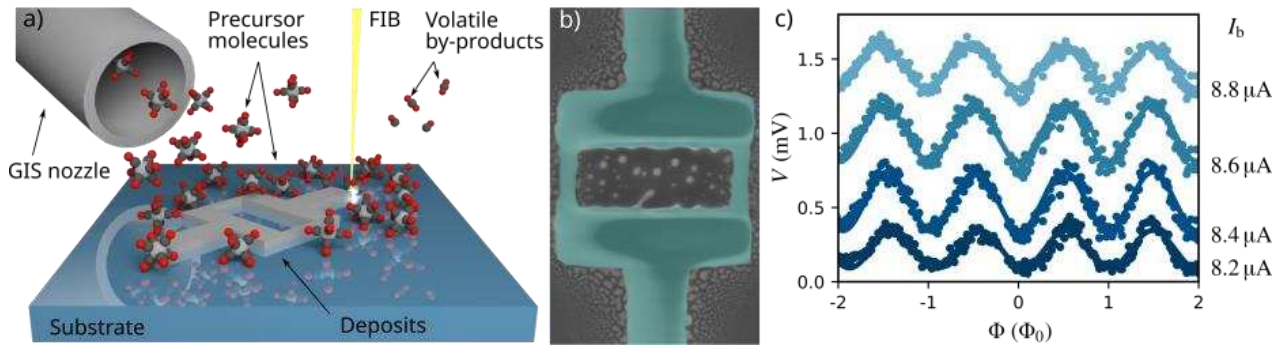


Fig. 1: a) Schematic representation of the FIB induced deposition of a nanostructure. b) A SEM image of a W-C (blue) nanoSQUID fabricated by  $\text{Ga}^+$  FIBID. c) The corresponding periodic modulation of the voltage,  $V$ , in dependence of the magnetic flux threading the SQUID loop,  $\Phi$ . The different lines correspond to different bias currents  $I_b \sim I_c$ . Reproduced from [3].

# Comparing Contrast in nano-CT and FIB/SEM Tomography of an Al Cast Alloy – A Correlative Microscopy Study

C. Pauly<sup>\*1</sup>, M. Engstler<sup>1</sup>, J. Fell<sup>2</sup>, H.-G. Herrmann<sup>2</sup>, F. Mücklich<sup>1</sup>

<sup>1</sup> Saarland University, Functional Materials, Campus D3 3, 66123 Saarbrücken, Germany

<sup>2</sup> Saarland University, Lightweight Systems, Campus E3 1, 66123 Saarbrücken, Germany

\* corresponding author email: [c.pauly@mx.uni-saarland.de](mailto:c.pauly@mx.uni-saarland.de)

Aluminum alloys with silicon as the main alloying element are widely used materials for general machine and automotive parts. The microstructure of hypereutectic Al-Si-alloys (>12 wt.% Si) comprises primary Si particles and a complex Al-Si eutectic. Elements such as Cu, Mg and Ni are added to form hard intermetallic precipitates which improve the mechanical properties at elevated temperatures. Size, shape and connectivity of these complex-shape intermetallic particles play a decisive role for the properties and are thus of special interest. As classic 2D methods cannot capture connectivity, a 3D analysis of the microstructure is necessary. Due to the particle size in the  $\mu\text{m}$  range, high-resolution 3D imaging techniques are required.

The x-ray computed tomography (XCT) has made strong progress in the past years enabling a resolution of  $\sim 50$  nm in lab-scale instruments. In the field of focused ion beam (FIB) serial sectioning, the introduction of the xenon plasma source has significantly increased the upper volume to an edge length of  $\sim 100\text{-}200$   $\mu\text{m}$ . Thus, both techniques overlap in terms of volume and resolution and correlative studies become possible.

In this work we use both XCT and FIB serial sectioning to image the microstructure of an AlSi13 cast alloy. Data was partly acquired from the same sample volume to compare the information given by the different techniques. Pros and cons of both techniques are discussed.

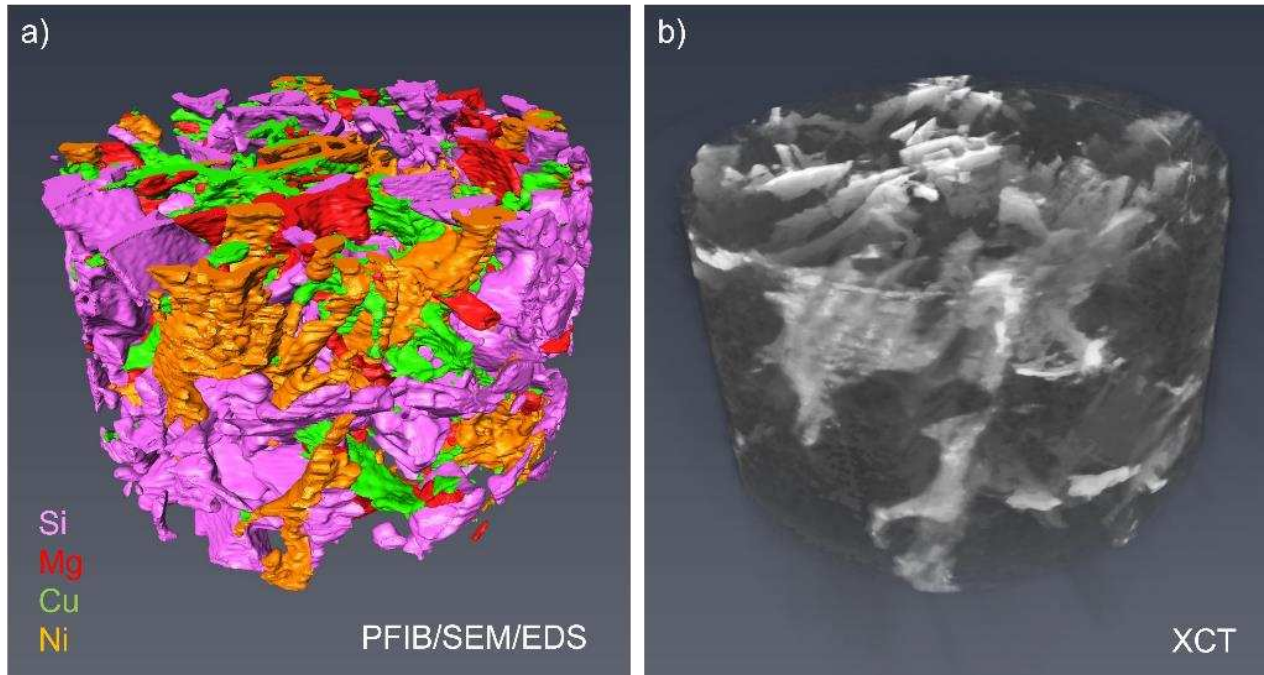


Figure 1: a) Reconstructed phases based on EDS data. Al matrix phase not shown for clarity, coloring according to principal elements. b) Volume rendering of XCT absorption contrast data. Phases containing heavy elements (Ni, Cu) appear bright.

# Explore unlimited process pathways for FIB nanopatterning and ion imaging using VELION

T. Richter, J. Stodolka\*, A. Nadzeyka

Raith GmbH, Konrad-Adenauer-Allee 8, 44263 Dortmund, Germany

[torsten.richter@raith.de](mailto:torsten.richter@raith.de)

Focused Ion Beam (FIB) nano patterning has become established as a versatile and precise fabrication method of manifold features at the nanoscale. Applications in nanoscale science require high resolution fabrication techniques at high fidelity, accuracy and reproducibility over multiple write fields in an automated manner.

VELION's configurable multi-ion species FIB technology enables tailoring of various nanostructures according to application related challenges. Various ion species can be selected from universal ion sources providing fast or slow and light or heavy ions from a single source [1]. This approach paves the way for unlimited process pathways based on a FIB system merged with a true lithography platform. Unmatched large-area FIB patterning and unlimited perfect write field stitching or patterning overlay facilitates various patterning strategies according to specific applications. As VELION utilizes comprehensive automation for unattended, uninterrupted reliable nanofabrication over several days, the instrument is the ideal companion for machine driven nanofabrication.

In this contribution we present the instrument concept, an overview of various nanofabrication approaches and application such as direct patterning, hard masking, or sample functionalization.

Beyond nanopatterning, VELION FIB with its Liquid Metal Alloy Ion Sources (LMAIS) provides excellent ion beam imaging capabilities [2]. Lithium is the lightest ion for LMAIS available from periodic table and provides sub 2nm lateral image resolution. Latest results of 3D Mill & Image sample analysis will be presented utilizing the best depth milling resolution with Bismuth and superior lateral resolution with Lithium ions.


[1] J. Gierak, P. Mazarov, L. Bruchhaus, R. Jede, L. Bischoff, Review of electrohydrodynamical ion sources and their applications to focused ion beam technology, JVSTB 36, 06J101 (2018).

[2] N. Klingner, G. Hlawacek, P. Mazarov, W. Pilz, F. Meyer, and L. Bischoff, Imaging and milling resolution of light ion beams from helium ion microscopy and FIBs driven by liquid metal alloy ion sources, Beilstein J. Nanotechnol. 11, 1742 (2020)




*Left: General VELION setup with vertical FIB column, Laser Interferometer Stage and SEM*


*Right: Multi ion species FIB from Liquid Metal Alloy Ion Sources.*

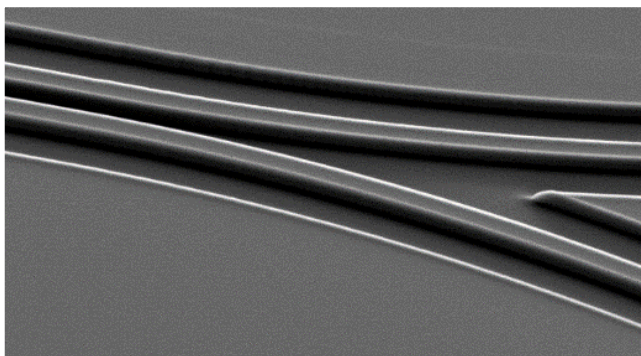


*Waveguide coupler*

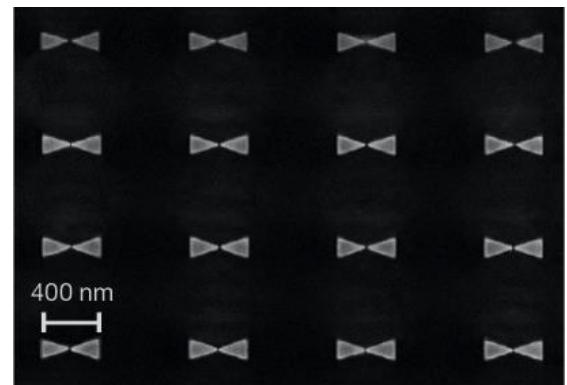

*Stepwise bowtie fabrication with Bi and Li ions*

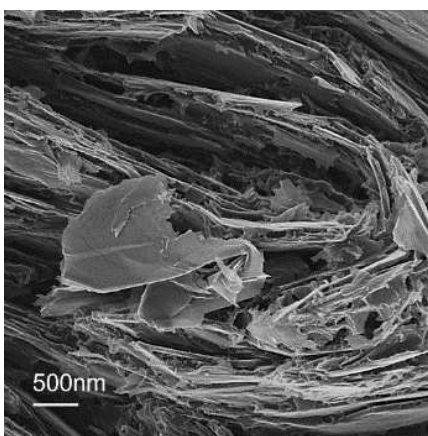

*Lithium-ion image of graphite*

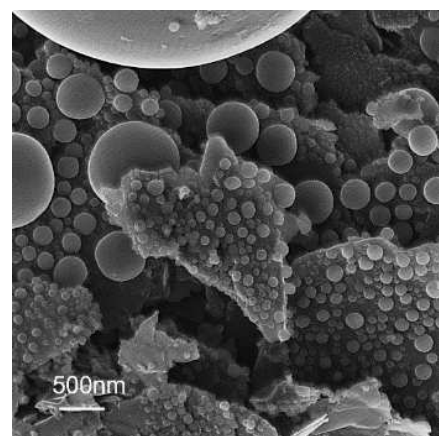

*Lithium-ion image of Sn/C*



# 3D Nanoprinting of Electrical AFM Nanoprobes

L. M. Seewald<sup>1</sup>, R. Winkler<sup>1</sup>, M. Brugger-Hatzl<sup>2</sup>, G. Kothleitner<sup>2,3</sup> and H. Plank<sup>1,2,3</sup>

<sup>1</sup> Christian Doppler Laboratory – DEFINE, Graz University of Technology, 8010 Graz, Austria

<sup>2</sup> Institute of Electron Microscopy and Nanoanalysis, Graz University of Technology, 8010 Graz, Austria

<sup>3</sup> Graz Center for Electron Microscopy, 8010 Graz

\* corresponding author email: [harald.plank@felmi-zfe.at](mailto:harald.plank@felmi-zfe.at)

Along other Scanning Probe Microscopy techniques Atomic Force Microscopy (AFM) has evolved into a widely used characterization technology at the micro- and nanoscale. Highest lateral resolution is enabled using very sharp probes which are scanned across a sample. This approach allows to surpass Abbe's limit of diffraction and has thereby paved the way towards new insights and innovative applications in research and development. In addition to high resolution topo-graphical imaging, AFM enables simultaneous mapping of surface properties such as mechanical, electrical, chemical, magnetic or thermal, with nanometer resolution. Advanced operation modes, however, rely on functionalized probes. Commonly, standard Si tips are coated with relevant materials to induce a desired functionality. Disadvantages are a larger apex radius, reducing resolution capabilities, and the risk of delamination effects due to mechanical stress during AFM operation, which reduce or even eliminate the targeted sensitivity. Hence fully functional uncoated probes with very sharp apex radii would be desired. With its unrivaled flexibility in terms of structural delicacy, geometrical flexibility and minimal requirements to substrate material and morphology, Focused Electron Beam Induced Deposition (FEBID)<sup>[1]</sup> is perfectly suited for the fabrication of AFM probes. Post-growth treatments enable accurate tailoring of material's properties towards the intended application<sup>[2]</sup>. With this motivation in mind, we here present a Pt-based 3D hollow cone concept for application in electrical AFM modes (CAFM, EFM, KPFM). The presentation will cover the entire interlinked design and fabrication process to achieve mechanically stable probes with sub-10 nm apex radii as well as their chemical transfer into highly crystalline Pt structures preserving the mentioned shape aspects (Fig.1b). We then present AFM studies to compare the performance of our FEBID probes to commercially available probes both in terms of resolution (Fig.1c,d) and functionality (Fig.1e). Together with the fact, that those concepts are meanwhile patented in collaboration with our industrial partners, this contribution clearly shows the industrial relevance of 3D-FEBID in the area of Atomic Force Microscopy.

[1] Winkler et al, *3D nanoprinting via focused electron beams*, J. Appl. Phys. 125, 210901 (2019)

[2] Geier et al., *Rapid and Highly Compact Purification*, J. Phys. Chem. C 118, 14009 (2014)



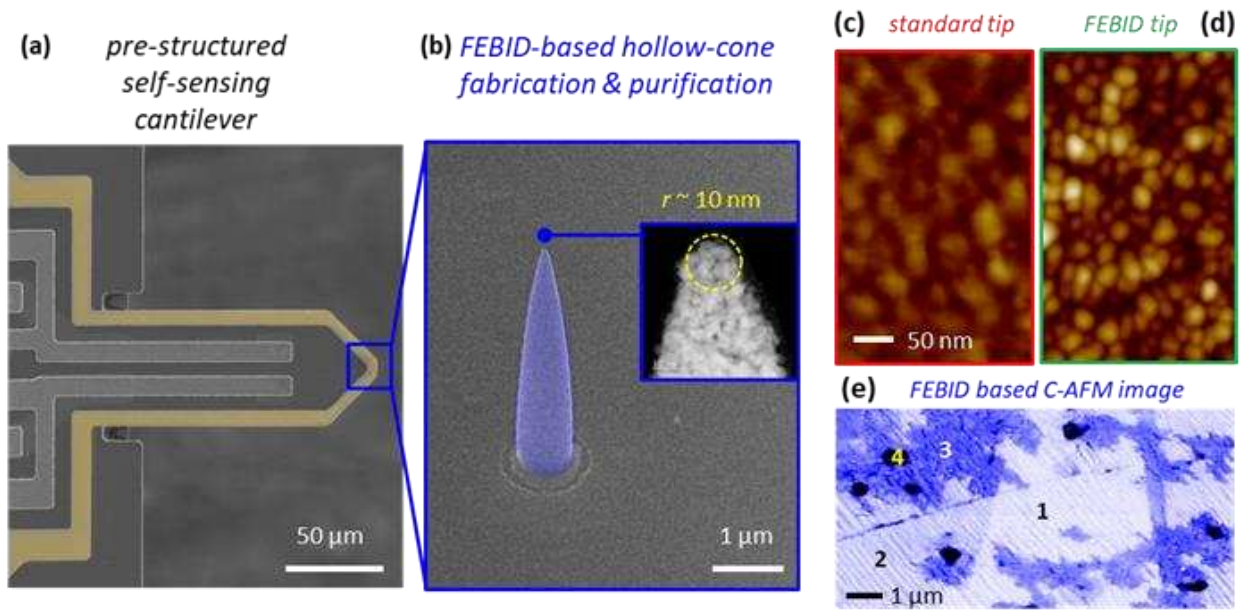


Figure 1: (a) pre-structured AFM self-sensing cantilever, which was first equipped with a PtC hollow-cone (b), further transferred into pure Pt with apex radii in the sub-10 nm regime (inset in b). (c) and (d) give a direct AFM topography comparison between commercially available C-AFM probes and the here relevant hollow-cones, respectively, which clearly shows the improved lateral resolution. (e) representative CAFM image enabling clear identification of a single-layer (1) and multi-layer graphene (2), copper oxide particles (3) and insulating Antimony particles (4) and proves CAFM applicability.

# **Manipulation and Study of Antiferromagnetic Order Enabled by Focused Ion Beam Fabrication**

S. C. Haley<sup>1,2</sup>, E. Maniv<sup>1,2,3</sup>, and J. G. Analytis<sup>1,2\*</sup>

<sup>1</sup> Department of Physics, University of California, Berkeley, California 94720, USA

<sup>2</sup> Materials Sciences Division, Lawrence Berkeley National Laboratory, Berkeley, California 94720, USA

<sup>3</sup> Department of Physics, Ben-Gurion University of the Negev, Beer-Sheva 84105, Israel

\* corresponding author email: [analytis@berkeley.edu](mailto:analytis@berkeley.edu)

In ferromagnetic solids, electron spins prefer to align with each other, with macroscopically observable results and straightforward applications: souvenirs can stick to refrigerators, stable magnetic fields can translate electric currents into sound, and information can be stored in the direction of the collective magnetic moment of a material. In antiferromagnets, on the other hand, spins prefer to be anti-aligned with their neighbors. In practice, this yields much more complicated and varied patterns on a microscopic level, which are more challenging to study directly because they do not generate a net external magnetic field and are difficult to manipulate with one. We are nonetheless interested in these systems, however, both for their exotic phase diagrams and low-energy excitations, and for their use in technological applications made possible by their complex magnetic properties.

By fabricating specialized transport devices using the FIB, we have explored and tested a first-generation prototype for low-power computing components based on antiferromagnets. These devices have also found value as a novel tool for identifying unusual magnetic dynamics and textures that are challenging to study by other means.

# **Backside FIB/SEM analysis strategy to identify a new failure mode at an automotive magnetic sensor device**

M. Simon-Najasek<sup>\*1</sup>, S. Hübner<sup>1</sup>, M. Lejoyeux<sup>1</sup>, A. Lindner<sup>2</sup>, F. Altmann<sup>1</sup>

<sup>1</sup> Fraunhofer Institute for Microstructure of Materials and Systems IMWS, 06120 Halle

<sup>2</sup> TDK-Micronas GmbH, Hans-Bunte Strasse 10, 79108 Freiburg

\* corresponding author email: [michel.simon-najasek@imws.fraunhofer.de](mailto:michel.simon-najasek@imws.fraunhofer.de)

Integrated circuits (IC) used in automotive applications must be produced under high quality standards. To achieve this goal each single fail from qualification or also field application must be traced down to its physical failure root cause to increase the overall device reliability. In most cases the failure analysis (FA) strategy of IC related defects is straight forward by electrical testing, followed by defect localization and finally destructive physical analysis to identify the defect structure and derive the root cause of the electrical malfunction.

In some cases, such a strategy is not straight forward, especially if the failure site can not be exactly located or if the IC-package interaction causes the failure mode.

In this work we will demonstrate a failure root cause analysis of an automotive magnetic field sensor IC which fails by an electrical short between 2 metal layers, fig. 1. Over time a conductive path is created and at a certain resistance the devices fail with increased power consumption. This could be shown several times, but the question stays unanswered over a longer time regarding the primarily root cause. There are different hypotheses discussed, for example electrostatic discharge, package issues as well as imprints induced by probe tools.

In a standard FA procedure, the package is fully removed above the IC to get access to its functional surface structure. Then further localization and preparation steps are following. At this specific case the original root cause could not be identified from frontside, that's why it was decided to turn the preparation starting from backside without removing the package at the frontside. The real challenge was that the failing position could not be localized in detail. Due to the electrical behavior of the sensor and the experiences collected by many fails before the position could be estimated and correlated to a several 10µm long metal network.

To get access to the ROI the Si substrate was mechanically thinned down to around 30µm thickness. After this the remaining silicon was removed by Xe-Plasma FIB trenching down to the active IC structure, fig. 2. Then the metal network was cut

by Ga-FIB at several positions and investigated by passive voltage contrast to identify the shorted part, fig. 3. Finally, the remaining shorted metal lines were screened by forwarding FIB cross sectioning under SEM observation till irregularities (delamination, cracks) could be found in the IC structure. The defect origin could be identified due to the preserved mould compound on top of the IC. A large filler particle touches the IC passivation at the position of the cracks. Due to the mechanical stress during the moulding process these cracks were initiated, later metal was migrating into these cracks forming an electrically conductive path between metal 1 and 2, fig. 4. The same failure mode could be verified at several sensor devices.

Fig. 1: SEM image of a FIB cross-section showing shorted metal lines due to a metal filled crack in-between these lines

Fig. 2: P-FIB backside trenching of the Si substrate overlaid with the IC layout of the ROI, before and after detailed Ga-FIB investigation

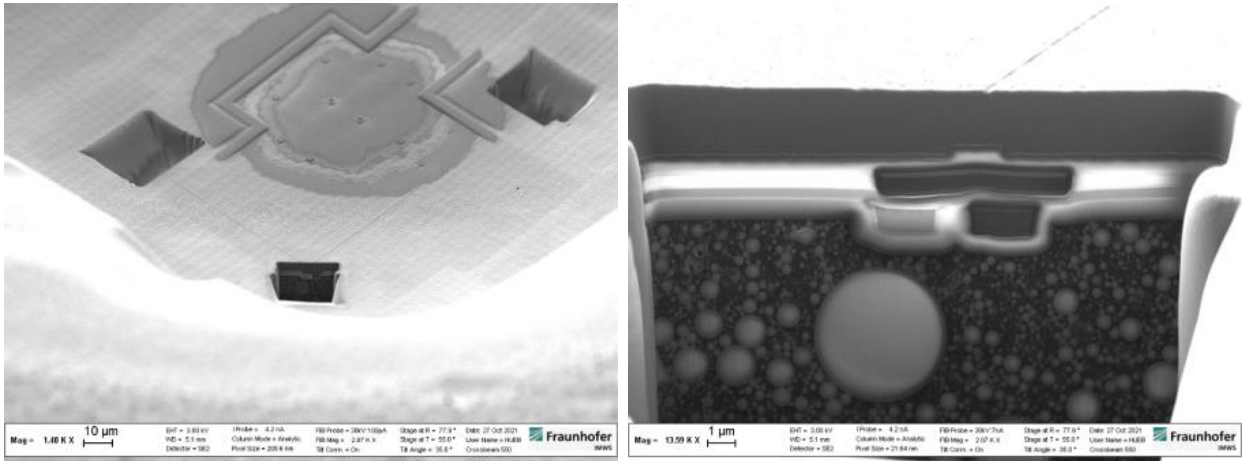


Fig.3: Passive voltage contrast investigation to isolate the defect position in the IC network

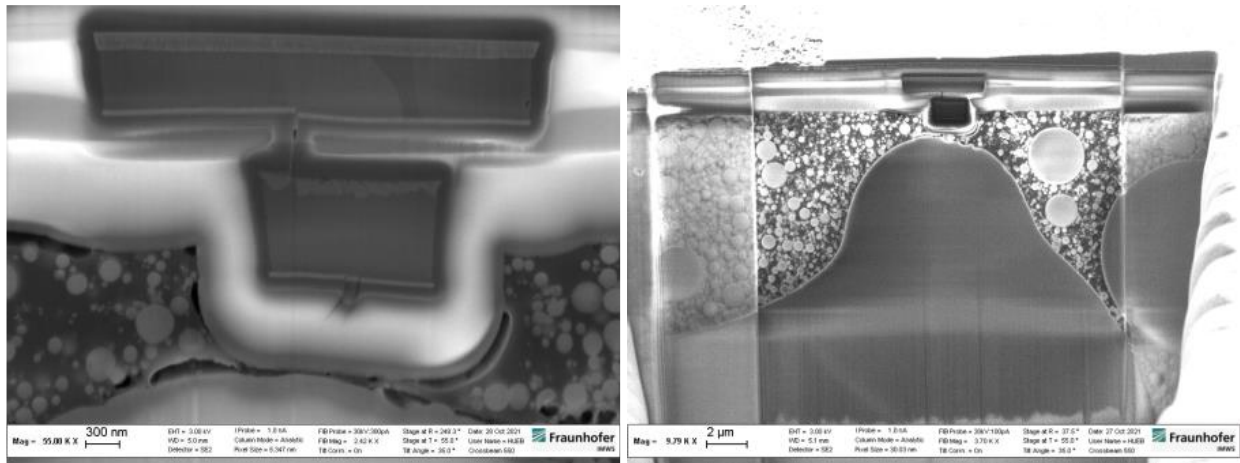


Fig.4: SEM cross-section showing a delaminated mould-passivation interface with cracked chip passivation and interlayer induced by a large filler particle

# **in situ Force Measurements - the FIB/SEM as a Mechanical Characterization Tool**

A. J. Smith\*, K. Schock, A. Rummel, S. Kleindiek

Kleindiek Nanotechnik, Aspenhastr. 25, 72770 Reutlingen, Germany

\* corresponding author email: [andrew.smith@kleindiek.com](mailto:andrew.smith@kleindiek.com)

Mechanically characterizing materials inside FIB/SEMs is a common task with a myriad of use cases in various fields of research.

In this work, the authors will present several examples of how in situ force measurements can be performed using different setups inside the FIB/SEM's chamber, thus utilizing the FIB/SEMs capability to modify, customize or otherwise prepare samples for testing as well as image samples from different angles for a more comprehensive set of images for later analysis.

Examples include bending FIB cut beams and comparing EBSD results obtained prior to and post bend (Fig. 1), flat punch experiments for elucidating forging properties of novel superalloys (Fig. 2), characterizing nanowires, CNTs, and other structures (Fig. 3), etc.

The described measurements are achieved using one of three setups: 1. Smallest forces - in the range of some nN - can be measured using self-sensing, piezo-resistive AFM cantilevers. 2. Another option is to utilize moveable sample holders with precisely calibrated spring constants. In this manner the deflection observed in the FIB/SEM can be used to calculate the applied force. By choosing from spring loaded sample holders with varying spring constants, a wider range of forces can be addressed. 3. Force transducers can be used to measure large forces up to several N.

Each approach has its own distinct use cases, advantages, and disadvantages. These will be discussed, as well.

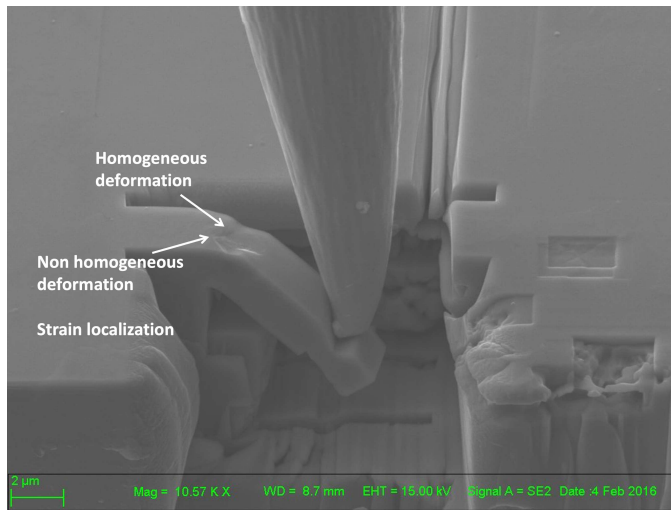


Fig. 1: FIB cut beam after bending failure (courtesy Archie, MPIE)

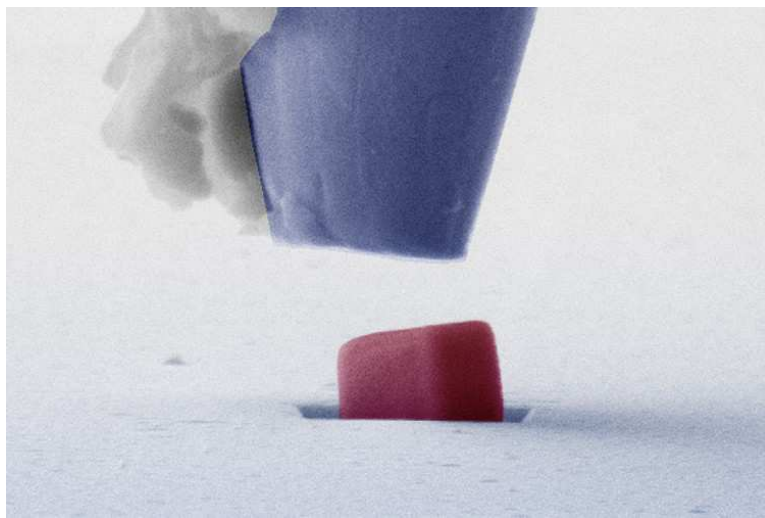


Fig. 2: Flat punch of a super alloy cube set in a FIB-cut base (courtesy Roesler, TU Braunschweig)

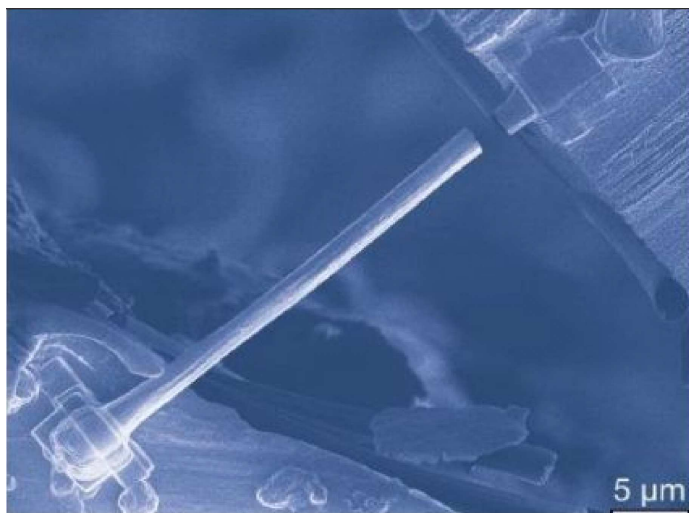


Fig. 3: Tensile test on a biological structure (courtesy Wegst, MPI Stuttgart)



# Development and SEM integration of the Nano Aperture Ion Source

M.L. Simons<sup>1\*</sup>, W.D. Laur<sup>1</sup>, T. van den Brink<sup>1</sup>, P. Kruit<sup>1</sup>, C.T.H. Heerkens<sup>1</sup>, A. Mahgoub<sup>1</sup>

<sup>1</sup> TU Delft


\* corresponding author email: [m.l.simons@tudelft.nl](mailto:m.l.simons@tudelft.nl)

The Nano-Aperture Ion Source (NAIS) is a new type of ion source, researched by L. van Kouwen, which should enable several ion species ( $\text{Ga}^+$ ,  $\text{He}^+$ ,  $\text{Ne}^+$ ,  $\text{Xe}^+$ ) without compromising brightness<sup>1</sup>  $\geq 1 \cdot 10^5 \text{A}/(\text{m}^2 \text{srV})$ , or energy spread. An electron beam from a Schottky source is focused on a hole in a nano gas channel. Gas flowing through this channel is ionized due to electron impact ionization, after which it is extracted from the source by an extraction potential (*Fig.1*).<sup>2</sup> L. van Kouwen identified problems with chip design and experimental setup. The production yield of the NAIS chips, created by MEMS technology, was low because two membranes needed a voltage difference while maintaining a 1-micron separation. Therefore, the chip design has been altered. These changes have considerably improved production yield and the maximum pressure the chip can withstand, up to 3 bar. This alteration does negatively affect the possibility to extract the ions from the chip, since an extraction field can only be applied outside the chip.


Instead of an electric field within the chip, a structure can be used to guide the gas towards the extractor, which can be created using Electron Beam Induced Deposition (EBID) (*Fig.2*). To optimize this structure, several simulations are performed in the Rarefied Gas Regime using Python and SPARTA. Together with an electric field and ion trace simulations in GPT, a theoretical brightness can be acquired and optimized by changing the structural design. An experimental setup is being used to test chip life time and perform brightness measurements (*Fig.3*). This setup is modified by applying a new extractor with a 6D-aligner directly onto the chip and by decoupling the chip position and the electrostatic components. However, if the NAIS chip is to be used for FIB applications, a large demagnification of the source is required to reduce the ion beam spot size. Therefore, an additional setup is being developed by combining a commercially available SEM and FIB columns in one device on the same beamline (NiCole and Tomahawk columns provided by Thermo Fisher Scientific). This setup is scheduled to be ready in 2023 and should produce a usable prototype.

[1] Leon van Kouwen, Pieter Kruit; Brightness measurements of the nano-aperture ion source; Journal of Vacuum Science & Technology B (2018).


[2] Leon van Kouwen; Introduction to focused ion beams, ion sources, and the nano-aperture ion source; Advances in Imaging and Electron Physics Including Proceedings CPO-10 (2019), Start Page 181.

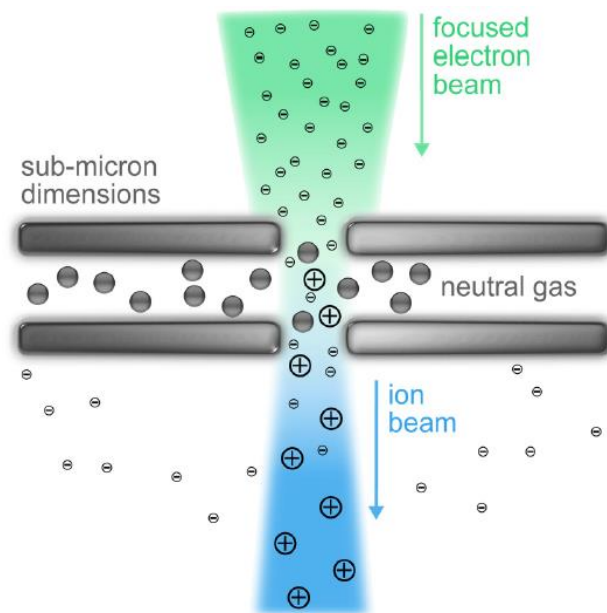

*Fig.1* principle of the NAIS: gas is ionized by a focused electron beam and extracted by a strong electric field.

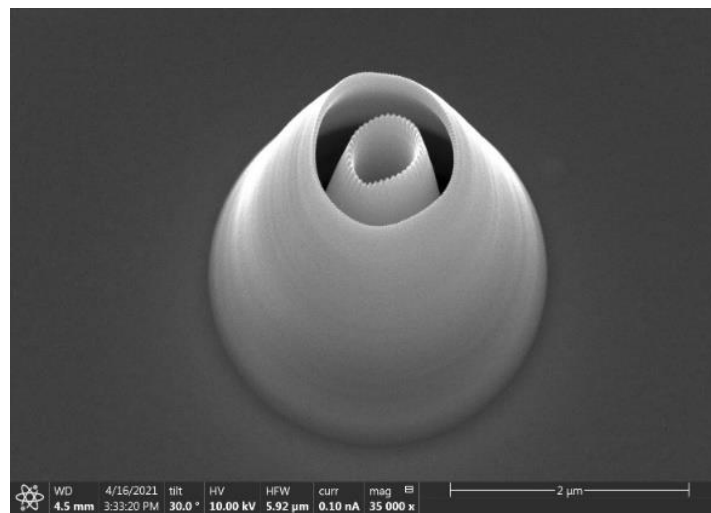

*Fig.2* produced microstructure (EBID) to enhance ion beam performance.

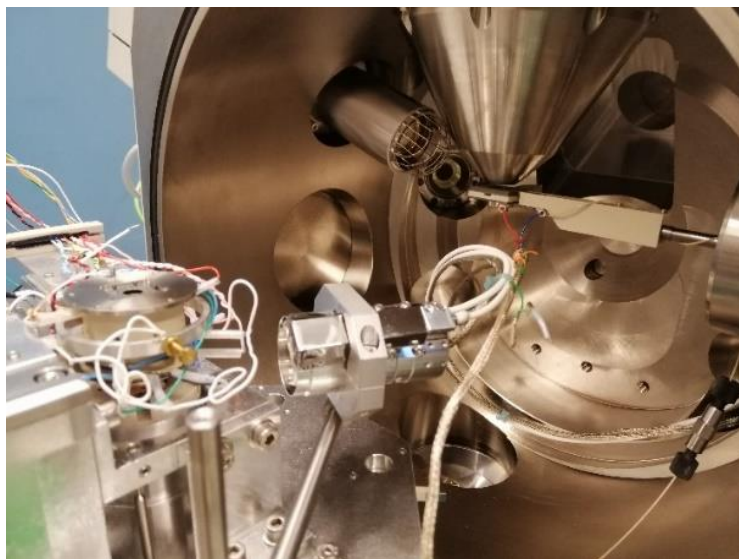

*Fig.3* current experimental setup in SEM chamber.

## **Extending microstructure characterization from nm to mm scale using the latest multiple ion species plasma FIB integrated with femtosecond laser**

Min Wu, Bartłomiej Winiarski and Remco Geurts

Thermo Fisher Scientific, Eindhoven, The Netherlands.

Conventionally 3D X-ray scanning technique such as microCT recovers large sample volume however suffers from relatively low voxel resolution, meanwhile high resolution 3D volume methods such as serial sectioning and imaging using a FIB-SEM can only recover 3D volume in the order of  $\leq 40 \times 40 \times 40 \mu\text{m}^3$ . Plasma FIB-SEM expands these techniques to volume  $\sim 250 \times 250 \times 250 \mu\text{m}^3$  keeping the voxels sizes in the dozens of nm-ranges. Recently, femtosecond Laser PFIB-SEM pushed these 3D techniques further to mm-scale volumes, setting the standards for multi-modal data collection from nm to mm scales and bridge the gap between microCT and FIB-SEM, while maintaining the advantage of high resolution SEM imaging of the cut face.

In this abstract we present the latest development of multiple ion species plasma FIB integrated with femtosecond laser, the Thermo Scientific Helios 5 Laser Hydra system. We address the effects of ultra-short pulse laser ablation and discuss various application usecases using this latest versatile technology. We will also discuss the possibility of laser ablation and slicing materials with multiple ion species under cryogenic conditions using our latest 360° rotatory Aquilos cryo stage, and the newly developed inert gas transfer system CleanConnect.

We will also present the latest automation software from Thermo Fisher Scientific, AutoTEM 5 for lamellae fabrication and Auto Slice & View 5 for 3D serial sectioning and volume reconstruction.

AutoTEM 5 facilitates users to prepare the highest quality S/TEM samples and cross-sectioning with a reliable fully-robotic fashion. It is a unique solution with highly configurable workflow to enable preparation of a wide range of samples for Inverted, Top Down and Planar use cases. The complete in-situ sample preparation workflow features fully automated, unattended multi-site chunk milling, lift-out and final thinning.

Auto Slice & View 5 is the latest generation of Thermo Fisher Scientific's automation software for automated serial sectioning and imaging through a user-defined volume of a specimen. Auto Slice & View 5 provides all necessary features for high precision slicing and imaging through reliable automation, easy setup, monitoring and on-the-fly correction of project parameters. Some of the notable new features include editable job templates to simplify and speed up job setting, rocking beam polishing, and fast and reliable auto functions (auto focus, auto stigmation and auto lens alignment). In addition, a newly developed Spin mill function allows the software to automate sequential milling of large (up to 1 mm) areas. The sample surface is oriented almost parallel to the ion beam and the stage is periodically rotated to a series of pre-defined positions until a full rotation of 360° is fulfilled. Multiple areas for image acquisition could also be configured per every Spin mill site.

# The patterning toolbox FIB-o-mat: Exploiting the full potential of focused ion beams for nanofabrication

V. Deinhart<sup>1,2\*</sup>, L.M. Kern<sup>3</sup>, J. Kirchhof<sup>4</sup>, S. Juergensen<sup>4</sup>, K.I. Bolotin<sup>4</sup>, S. Reich<sup>4</sup>, B. Pfau<sup>3</sup> and K. Höflich<sup>1,2</sup>

<sup>1</sup> Ferdinand-Braun-Institut gGmbH, Leibniz-Institut für Höchstfrequenztechnik, Gustav-Kirchhoff-Str. 4, 12489 Berlin, Germany

<sup>2</sup> Corelab Correlative Microscopy and Spectroscopy, Helmholtz-Zentrum Berlin für Materialien und Energie, Hahn-Meitner-Platz 1, 14109 Berlin, Germany

<sup>3</sup> Max Born Institute for Nonlinear Optics and Short Pulse Spectroscopy, Max-Born-Straße 2A, 12489 Berlin, Germany

<sup>4</sup> Freie Universität Berlin, Arnimallee 14, 14195 Berlin, Germany

\* corresponding author email: [victor.deinhart@fbh-berlin.de](mailto:victor.deinhart@fbh-berlin.de)

Future breakthroughs in nanotechnology will rely on the ability to fabricate materials and devices by design, i.e. to tailor both material properties and device geometries according to a sophisticated blueprint. Direct writing using focused beams of ions or electrons is a powerful technique not only for rapid-prototyping of novel device components but also for mask-less processing of delicate nanostructures and local modification of materials. To achieve an optimal result in the patterning process, a full control over the beam path including its rasterization is necessary. Here, we present the open-source Python toolbox FIB-o-mat [1] for automated pattern creation and optimization (cf. Fig. 1).

Patterning with an ion beam is a digital process in which the beam spot dwells for a defined time at a certain location and is then shifted by a defined distance ('pitch') to dwell again. This beam path can be created from a high-level geometry that is rasterized following a pre-defined routine. Therefore, FIB-o-mat provides various geometries including their Boolean combinations together with all possible rectangular and annular rasterization schemes. For complex geometries however, rasterization may lead to artifacts in narrow regions or for large curvatures. Here, FIB-o-mat offers a low-level approach where the beam path is generated by curve off-setting for any given input geometry to achieve the best patterning fidelity.

The functionalities of the Python-based toolbox FIB-o-mat are showcased for He ion beam processing of three different material systems. The magnetic properties of Co-based multilayers were locally modified using the high-level beam path generation combined with automation via stage control [2]. Plasmonic tetramer antennas were cut from single-crystalline gold on glass, demonstrating the ultimate patterning resolution of the focused He ion beam with gap sizes down to 3 nm (cf. Fig. 2). Apart from the low-level beam path generation, a local dose optimization had to be carried out. Finally, suspended single layer graphene was patterned into simple trampoline resonators but also complex phononic crystal structures [3] using the low-level beam path generation and automation via stage control.

- [1] V. Deinhart, L.M. Kern, J.N. Kirchhof, S. Juergensen, J. Sturm, E. Krauss, Th. Feichtner, S. Kovalchuk, M. Schneider, D. Engel, B. Pfau, B. Hecht, K.I. Bolotin, S. Reich and K. Höflich; *The patterning toolbox FIB-o-mat: Exploiting the full potential of focused helium ions for nanofabrication*; Beilstein J. Nanotechnol. 12 (2021), 304.
- [2] L.M. Kern, et al.; *Deterministic Generation and Guided Motion of Magnetic Skyrmions by Focused He<sup>+</sup>-Ion Irradiation*, Nano Lett. 22 (2022), 4028.
- [3] J. N. Kirchhof, et al.; *Tunable Graphene Phononic Crystal*; Nano Lett. 21 (2021), 2174.

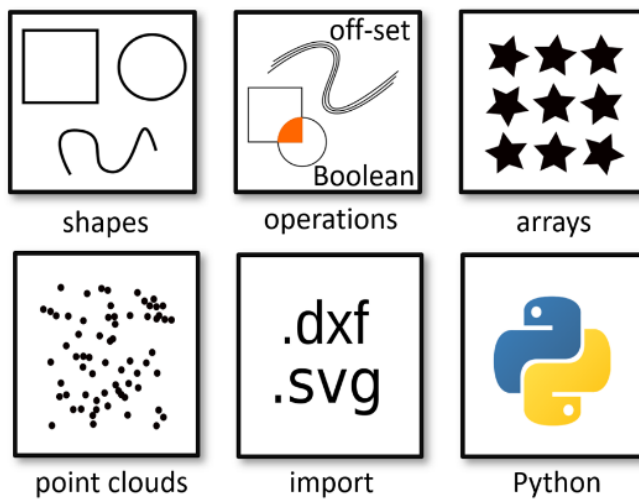


Fig. 1: FIB-o-mat tools

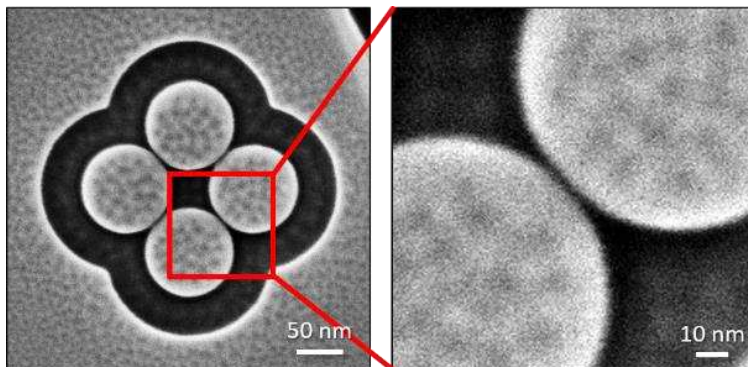


Fig. 2: Plasmonic tetramer antenna with ultimate gap resolution achieved by He ion beam patterning

# Focused Ion Beam Induced Strain Generation in Silicon Membranes

D. C. Cox<sup>\*</sup>, and M.G. Masteghin

National Ion Beam Centre, University of Surrey, Guildford Surrey, GU2 7XH, U.K.

<sup>\*</sup> corresponding author email [d.cox@surrey.ac.uk](mailto:d.cox@surrey.ac.uk)

In this presentation we show an incredibly simple and controllable way of generating strain in suspended Si membranes using a Focused Ion Beam. The strain we observe approaches record levels observed in this material while retaining high crystallinity[1].

What this study also shows however is the power of serendipity. Without doubt the work contained here would not have occurred were it not for very substantial luck. The FIB induced strained region does not occur in the ion irradiated regions but in the non-irradiated regions adjacent to it. It was only possible to observe this effect due to presence of existing compressive stress, leading to buckling, in our 35 nm thick membrane substrates. This combined with an experiment that required us to irradiate them with a 30 kV ion beam led to the phenomena reported here.

FIB Irradiation of samples with energetic ion beams can create many effects such as doping, sputtering, chemical modification and damage creation. These effects depend greatly on ion energy, dose and species. The creation of damage and/or amorphisation due to the implantation of energetic ions is well understood and none more so than in the silicon system due to its prominence in the semiconductor industry. For a normal incidence Ga ion at 30 kV the average ion range can be calculated to be approximately 25 nm deep, or a little less for an equivalent energy Xe ion. In both cases full amorphisation of the Si to a depth of the ion range occurs at less than  $4 \times 10^{14}$  ions/cm<sup>2</sup>. For a focused ion beam to achieve doses of this intensity is trivial taking only a few seconds even at modest beam currents.

If 20 kV Ga ions or 30 kV Xe ions are implanted into a 35 nm thick Si membrane the ion range and hence crystal amorphisation corresponds to very close to the midpoint of the membrane thickness. This creates a bi-layer system that induces bending of the implanted membrane regions straining the adjacent unimplanted regions. Using simple FIB patterning strategies and exploiting crystal anisotropy we have shown it is possible to create biaxial strain in excess of 3% and uniaxial in excess of 8%. If the same technique were applied to germanium a transition to a direct band-gap semiconductor could be achieved at equivalent strains.

**[1] Stress-strain engineering of single-crystalline silicon membranes by ion implantation: Towards direct-gap group-IV semiconductors**

M.G. Masteghin, V. Tong, E.B. Schneider, C.C. Underwood, T. Peach, B. N. Murdin, R.P. Webb, S.K. Clowes and D.C. Cox. *Physical Review Materials*, Vol.5, 124603. 2021

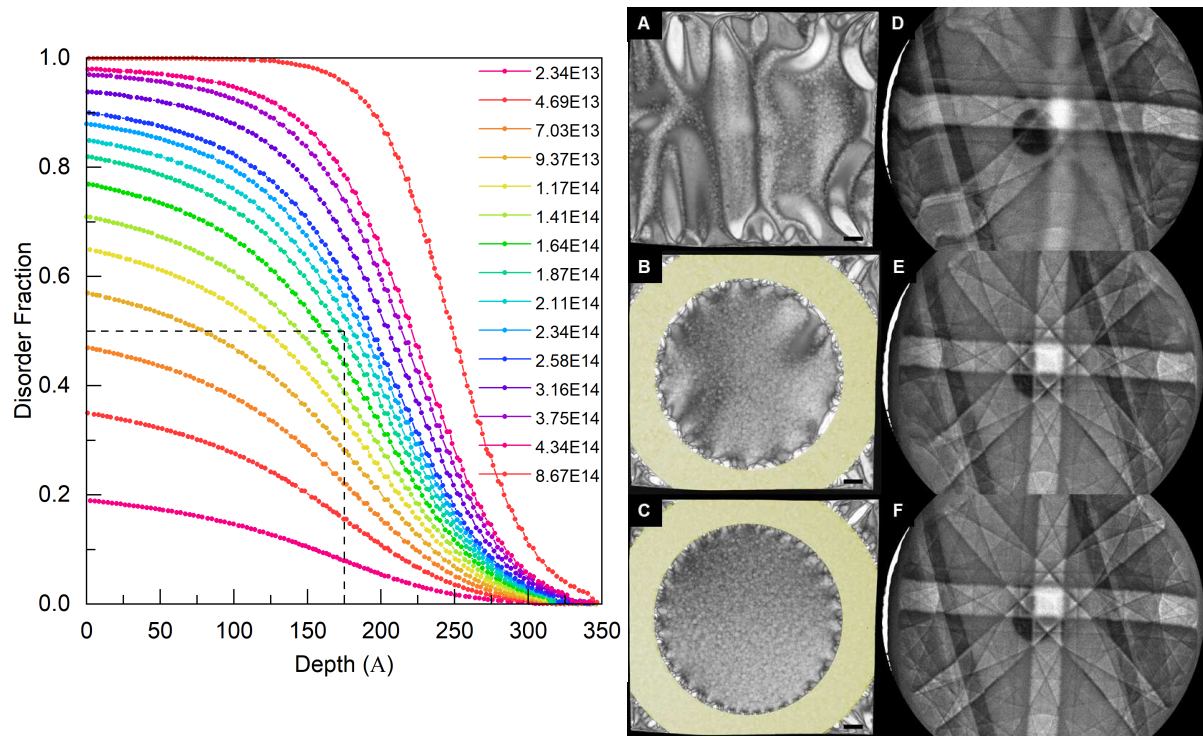


Fig. 1: (Left hand panel) Disorder fraction vs depth for 30 kV Xe ions into Si for the dose range  $2.4 \times 10^{13}$  to  $8.7 \times 10^{14}$  ions/cm<sup>2</sup> showing how the irradiated side of the membrane becomes amorphous but the lower side remains undamaged. (Right hand panel) The Si membrane prior to implantation (A), after partial exposure of the annular pattern, false coloured in yellow (B) and the same sample at a higher dose where the membrane buckles have been removed (C). The diffraction patterns demonstrate how the central region of the membrane is pulled flat but retains its crystallinity. The inner annular diameter is 80 µm.



# Two Microscopes are better than One – *In-situ* Correlative Analysis by combination of AFM, SEM, and FIB

H. Frerichs<sup>1\*</sup>, S. Seibert<sup>1</sup>, L. Stühn<sup>1</sup>, M. Wolff<sup>1</sup> and C.H. Schwalb<sup>1</sup>

<sup>1</sup> Quantum Design Microscopy GmbH, 64293 Darmstadt, Germany

\* corresponding author email: [frerichs@qd-microscopy.com](mailto:frerichs@qd-microscopy.com)

Combining different analytical methods in one instrument is of great importance for the simultaneous acquisition of complementary information. Especially the *in-situ* combination of scanning electron microscopy (SEM) and atomic force microscopy (AFM) enables completely new insights into the micro and nano-world. In this work, we present the unique *in-situ* combination of scanning electron / ion microscopy (SEM/FIB) and atomic force microscopy (AFM) for nanoscale characterization [1-2].

A particularly important aspect for advanced AFM measurements are the cantilevers itself. Especially the cantilever tips have to provide special properties to measure not only the topography, but also magnetic and electrical characteristics of materials. In this context, the FEBID process is a very promising approach for the design and fabrication of sophisticated tips with unique properties. We will present some case studies to highlight the advantages of interactive correlative *in-situ* nanoscale characterization for different materials and nanostructures, using FEBID-constructed tips. We show results for the *in-situ* electrical characterization by conductive AFM for 2D materials as well as electrostatic force microscopy (EFM) measurements of piezoceramic films that enables the precise analysis of grain boundary potential barriers in semiconducting BaTiO<sub>3</sub>-based ceramics [3]. The grain boundaries were located via BSE-SEM and measured afterwards using the *in-situ* EFM method. The barriers were shown to be significantly thinner and more pronounced as the amount of SiO<sub>2</sub> was increased from 0 to 5 mol% (see Fig. 1). These results can be directly correlated with electron backscatter diffraction (EBSD) measurements in order to link the AFM and SEM data to the crystallographic microstructure.

In addition, we will present results for the *in-situ* characterization of magnetic nanostructures by combination of SEM and high-vacuum magnetic force microscopy (MFM). For the *in-situ* MFM measurements, special high aspect ratio magnetic cantilever probes fabricated by electron beam induced deposition (FEBID) were used, which surpass conventional cantilevers in terms of lateral and magnetic resolution. The SEM enables to identify the grain boundaries on multilayer thin-film samples or stainless steel in order to measure the magnetic properties directly via MFM with nanometer resolution (see Fig. 2).

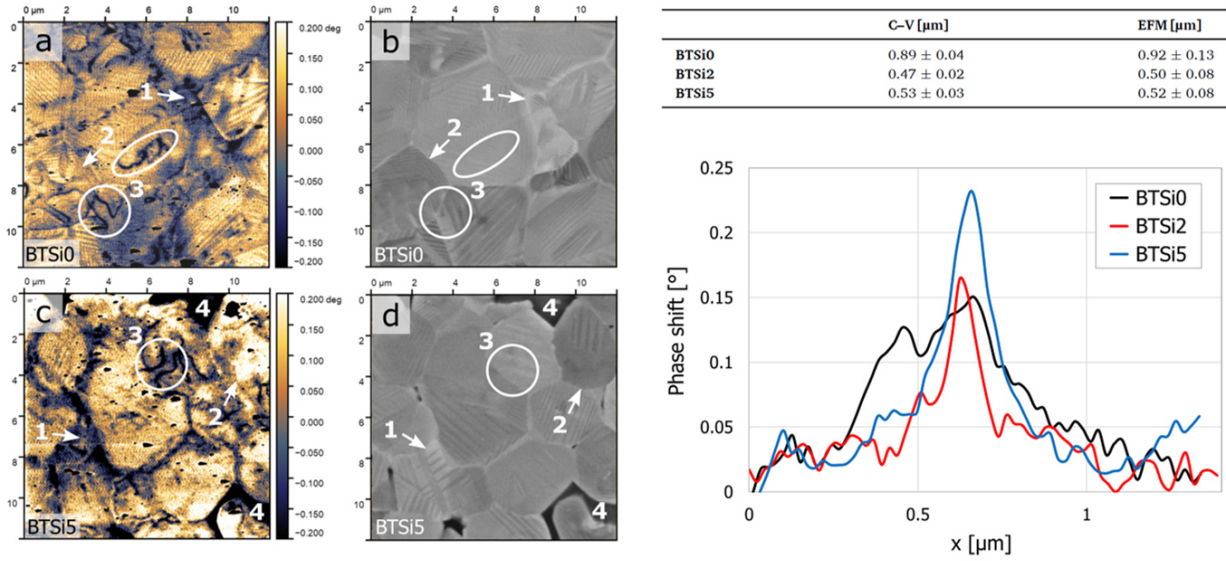


Fig 1: (Left) EFM signal of (a) BTSi0 and (c) BTSi5 with corresponding BSE-SEM images (b) and (d) of the exact same sample area. The EFM images were recorded with an external voltage of -3 V. (Right) Grain boundary potential barriers of BaTiO<sub>3</sub> based ceramics with varying SiO<sub>2</sub> content. The Images were extracted from [3].

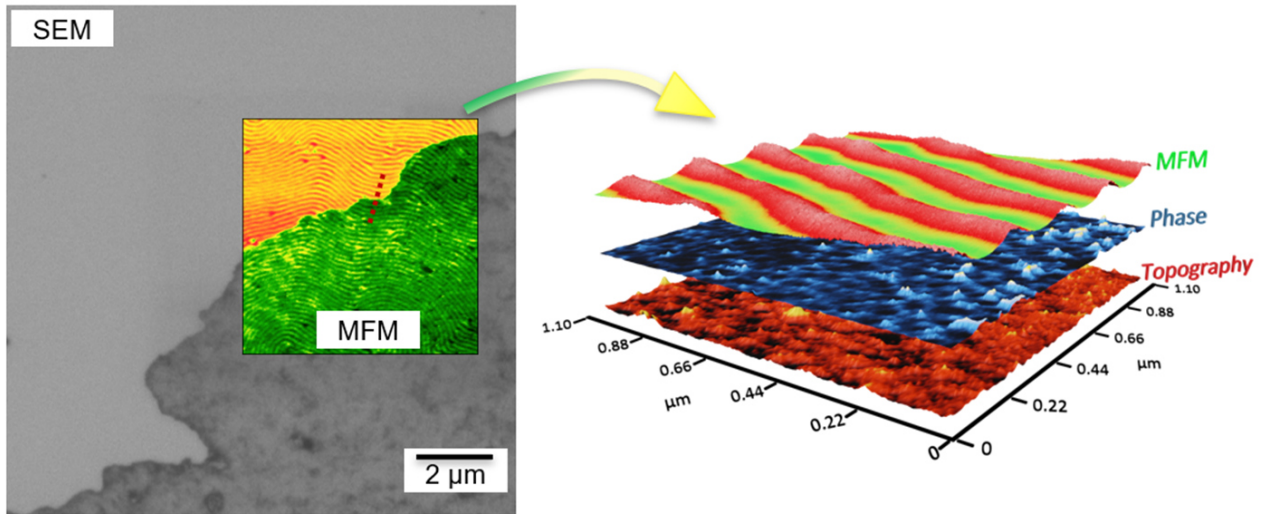


Fig 2: (Left) SEM image of a multilayer magnetic sample with correlative MFM image at a grain boundary. (Right) AFM data overlay of topography, phase, and MFM signals.

- [1] D. Yablon, et al.; *Cross-platform integration of AFM with SEM: Offering the best of both worlds*; Microscopy and Analysis 31 (2) (2017), 14.
- [2] S.H. Andany, et al.; *An atomic force microscope integrated with a helium ion microscope for correlative nanoscale characterization*; Beilstein Journal of Nanotechnology 11 (2020), 1272.
- [3] J.M. Prohning, J. Hütner, K. Reichmann, S. Bigl; *Complementary evaluation of potential barriers in semiconducting barium titanate by electrostatic force microscopy and capacitance-voltage measurements*; Scripta Materialia 214 (2022), 114646.

# Direct-Write 3D Nanoprinting of High-Resolution Magnetic Force Microscopy Nanoprobes

M. Brugger-Hatzl<sup>1</sup>, R. Winkler<sup>1</sup>, L. Seewald<sup>1</sup>, S. Mietsche<sup>2</sup>, T. Radlinger<sup>3</sup>, M. Huth<sup>4</sup>, S. Barth<sup>4</sup>, and H. Plank<sup>1,2,3,\*</sup>

<sup>1</sup> Christian Doppler Laboratory - DEFINE, Graz University of Technology, Steyrergasse 17, 8010 Graz, Austria

<sup>2</sup> Institute of Electron Microscopy and Nanoanalysis, Graz University of Technology, Steyrergasse 17, 8010 Graz, Austria

<sup>3</sup> Graz Centre for Electron Microscopy, Steyrergasse 17, 8010 Graz, Austria

<sup>4</sup> Institute of Physics, Goethe University, Max-von-Laue-Strasse 1, 60438 Frankfurt am Main, Germany

\* email: [harald.plank@felmi-zfe.at](mailto:harald.plank@felmi-zfe.at)

Magnetic devices play an important role in modern electronic, sensing or data storage applications. To exploit their full potential, high-resolution Magnetic Force Microscopy (MFM) is established as standard characterization technology as part of the research and development loop. Due to the ongoing trend towards smaller and smaller active feature sizes, the demands on high-resolution MFM tips are also increasing. Based on that motivation, we here aim on the fabrication of MFM nanoprobes with functional apex radii in the sub-10 nm regime. Traditional products mostly base on additional magnetic coatings, which increases the apex radii and therefore limits the lateral resolution during Atomic Force Microscopy (AFM) based MFM measurements. Another disadvantage of a magnetic coating is local delamination, which can occur due to the mechanical stress during scanning and lead to a change (or even complete loss) in magnetic sensitivity. Therefore, it was the goal to fabricate fully magnetic nanoscale tips, that do not require additional coating. Focused Electron Beam Induced Deposition (FEBID) was used for additive, direct-write 3D-nanoprinting of such magnetic tips on pre-finished self-sensing AFM cantilevers.<sup>[1],[2]</sup> For that, a novel  $\text{HCo}_3\text{Fe}(\text{CO})_{12}$  precursor was used, which is one of the few precursors, providing metal contents above 90 at.% after initial FEBID fabrication.<sup>[3]</sup> To explore the possibilities, we comprehensively studied the parameter space and their implications on morphology, structure and chemistry in detail by using SEM, EDX, and TEM and STEM EELS (Figure 1.a). Next, the tip geometry was further optimized by an advanced, dynamic patterning sequence to fulfil the high demands for AFM operation.<sup>[4]</sup> Additionally, the fabricated tips were subjected to different post-processing procedures such as post-irradiation with electrons, thermal treatments and purification protocols to explore and identify the most promising fabrication window. The basic performance of such MFM tips is then demonstrated with special focus on lateral resolution, magnetic phase shift and signal-to-noise ratio. Fully optimized FEBID-MFM tips were then tested on various magnetic samples (magnetic multilayer system (Figure 1.b-e), hard disc drives, magnetic recording tapes) and benchmarked to commercially available MFM tips (Figure 1.b-c). Finally, the wear resistance of such MFM nanoprobes was evaluated during a continuous operation scan over a period of 3.7 hours, which revealed the high durability of the presented concept (Figure 1.d-e). By that, we demonstrate the successful 3D-nanoprinting of MFM tips on self-sensing cantilevers, which fulfils the high requirements when aiming on industrially relevant MFM tips using FEBID-based 3D nanoprinting.

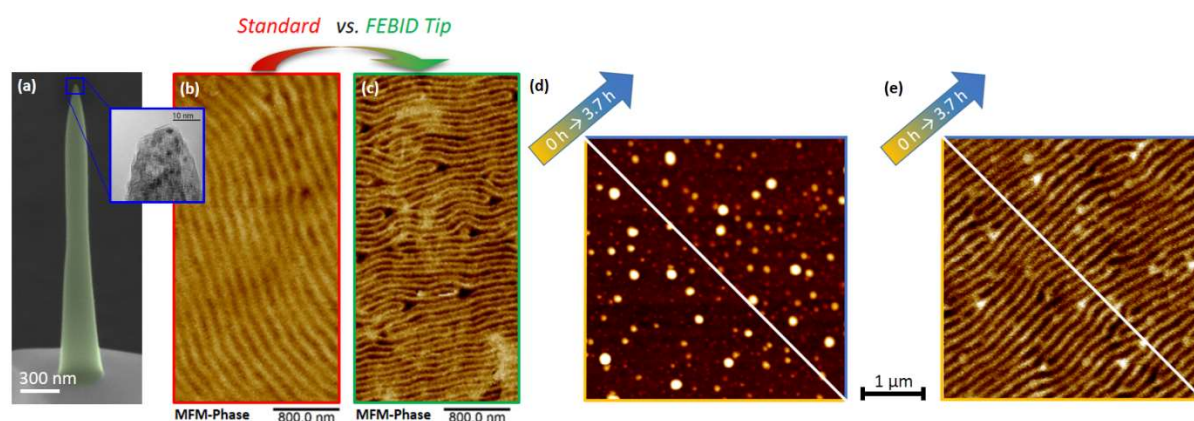


Figure 1: (a) SEM image of a specifically tailored,  $\text{Co}_3\text{Fe}$  3D nanoprobe for Magnetic Force Microscopy (MFM). The inset shows a TEM image of the tip region, revealing a fully crystalline tip radius around 10 nm. (b) and (c) show a direct comparison of magnetic MFM maps, taken with a commercial and a FEBID-based MFM tip, respectively, which reveal the superior performance of the latter. (d) and (e) demonstrate the wear resistance of FEBID tips via topography and magnetic phase maps, respectively. The lower left parts (yellow) show the starting situation, while the upper right parts (blue) give the results after continuous AFM / MFM operation of 3.7 hours, revealing practically identical results, which underlines the durability of FEBID based MFM tips.

## References

- [1] H. Plank et al., "Focused electron beam-based 3D nanoprining for scanning probe microscopy: A review," *Micromachines*, vol. 11, no. 1. MDPI AG, Jan. 01, 2020. doi: 10.3390/mi11010048.
- [2] J. Pablo-Navarro et al., "Magnetic Functionalization of Scanning Probes by Focused Electron Beam Induced Deposition Technology," *Magnetochemistry*, vol. 7, no. 10, 2021, doi: 10.3390/magnetochemistry7100140.
- [3] F. Porrati et al., "Direct writing of CoFe alloy nanostructures by focused electron beam induced deposition from a heteronuclear precursor," *Nanotechnology*, vol. 26, no. 47, p. 475701, Nov. 2015, doi: 10.1088/0957-4484/26/47/475701.
- [4] D. Kuhness et al., "High-Fidelity 3D Nanoprining of Plasmonic Gold Nanoantennas ", *ACS Appl. Mat. and Int.* 2021, 13, 1178-1191, doi: 10.1021/acsami.0c17030.

# Using the Helium Ion Microscope for Imaging and Modification of Nanostructures, 2D Materials, and SARS-CoV-2 infected Cells

Armin Götzhäuser

Physics of Supramolecular Systems and Surfaces,  
Bielefeld University, Universitätsstr. 25, 33615 Bielefeld, Germany

*e-Mail: ag@uni-bielefeld.de*

The Helium Ion Microscope (HIM) utilizes a focused beam of helium ions to image and modify materials with high spatial resolution, large depth of field, and chemical sensitivity [1]. HIM images show stronger chemical and topographical contrasts than images from the related scanning electron microscope, and the HIM is capable to resolve sub-nanometer features. Due to its charge compensation capability, the HIM can image insulating biological samples without additional conductive coatings. My presentation will contain examples of HIM imaging of nanomaterials, like 1 nm thick carbon nanomembranes (CNMs), 2D materials, and biological cells [2]. In an exploratory HIM study of SARS-CoV-2 infected Vero E6 cells, interactions between cells and virus particles, as well as among virus particles, could be imaged [3]. The HIM pictures show the three-dimensional appearance of SARS-CoV-2 and the surface of Vero E6 cells at a multiplicity of infection of approximately 1 with great morphological detail. The absence of a conductive coating allows a distinction between virus particles bound to the cell membrane and virus particles lying on top of the membrane. When applying higher ion currents, the HIM can be also used for the modification of materials. The capability of the HIM for nanolithography will be shown by milling 2D materials, where nanopores with diameters down to 1.3 nm were fabricated [4].

- [1] G. Hlawacek and A. Götzhäuser (Ed.): *Helium Ion Microscopy*, Springer-International (2016).
- [2] M. Schürmann, N. Frese, A. Beyer, P. Heimann, D. Widera, V. Mönkemöller, T. Huser, B. Kaltschmidt, C. Kaltschmidt and A. Götzhäuser: *Helium Ion Microscopy Visualizes Lipid Nanodomains in Mammalian Cells*, *Small* **43**, 5781 (2015).
- [3] N. Frese, P. Schmerer, M. Wortmann, M. Schürmann, M. König, M. Westphal, F. Weber, H. Sudhoff and A. Götzhäuser: *Imaging of SARS-CoV-2 infected Vero E6 Cells by Helium Ion Microscopy*, *Beilstein J. Nanotechnol.* **12**, 172 (2021).
- [4] D. Emmrich, A. Beyer, A. Nadzeyka, S. Bauerdick, J. C. Meyer, J. Kotakoski and A. Götzhäuser: *Nanopore Fabrication and Characterization by Helium Ion Microscopy*, *Appl. Phys. Lett.* **108**, 16310 (2016).

# Sample preparation and analysis of LLZO ceramics for solid state batteries with Cryo FIB/SEM and aberration corrected analytical STEM

C. Warres<sup>1</sup>, T. Lutz<sup>1</sup> and CJ. Burkhardt<sup>1\*</sup>

<sup>1</sup> NMI Natural and Medical Sciences Institute

\* corresponding author email: [claus.burkhardt@nmi.de](mailto:claus.burkhardt@nmi.de)

Solid state batteries are promising devices for high capacity energy storage in future electric vehicles. For stable, efficient and long life operation of these batteries stable interfaces of different layers are crucial. To optimize the interface between lithium and the solid electrolyte various interlayer materials are under development. Here we report our cryo-FIB/SEM, HRSTEM and EELS analysis of aluminum doped LLZTO (Lithium Lanthanum Tantalum Zircon Oxide) solid-state electrolyte coated with tin interlayer.

To analyze the morphology and chemical composition of coated LLZTO, FIB (focused ion beam) cross sections were prepared (Zeiss crossbeam 550 scanning electron microscope) and analyzed with EDS (Oxford X-Max 150). In a second step, a several  $\mu\text{m}$  thick lamella was prepared by FIB lift out technique [1]. To reduce beam induced damage the sample was finally polished inside the FIB/SEM under cryogenic conditions ( $-160^\circ\text{C}$ , Quorum PP3010). After finishing a TEM transparent lamella, the sample was warmed up to room temperature and transferred via an argon glove box on a TEM vacuum holder. High resolution analysis of the lamella was performed in a probe corrected JEOL ARM 200F electron microscope equipped with a JEOL Dual EDS system. The Li K line is not accessible in conventional EDS, therefor we used EELS (Electron Energy Loss Spectroscopy, Gatan Quantum ER Spectrometer) to measure the distribution of Lithium.

Results of STEM imaging and EELS mapping are shown in figure 1. STEM imaging with atomic resolution of the aluminum precipitate is shown in figure 2.

[1] Warres C., Lutz T. and Burkhardt C., Insight into LLZTO solid state battery material

[2] Mayo, M. and A.J. Morris, Structure Prediction of Li-Sn and Li-Sb Intermetallics for Lithium-Ion Batteries Anodes. Chemistry of Materials, 2017. 29 (14): p. 5787-5795



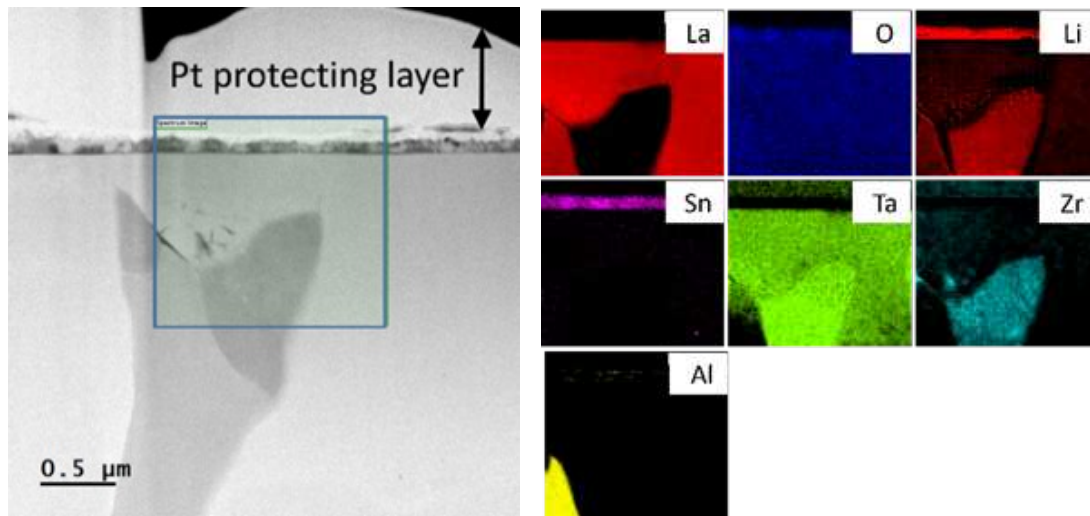


Fig. 1: STEM ADF image and EELS elemental mapping at 200keV in JEOL ARM200F.

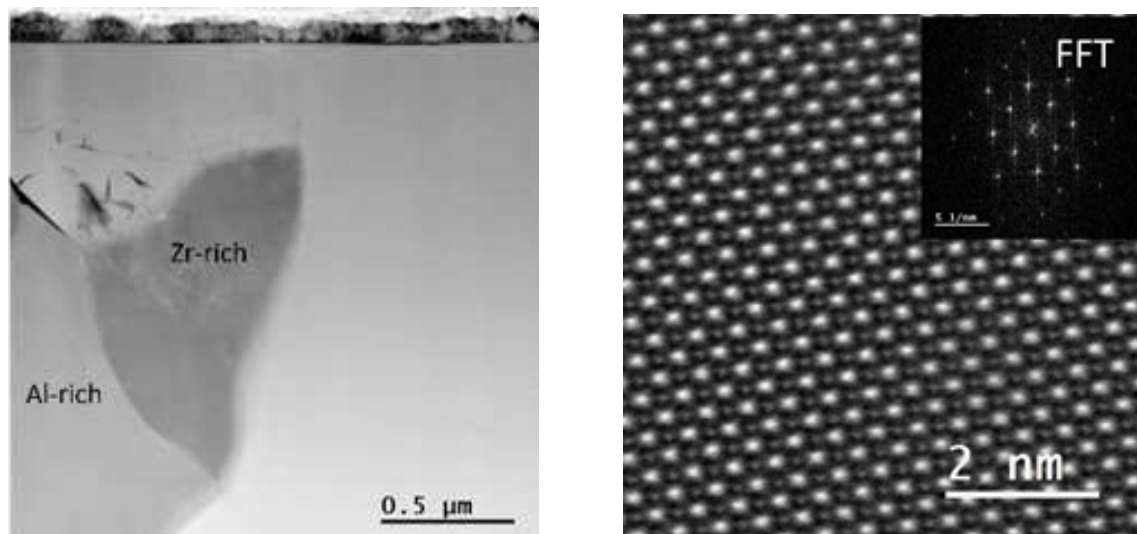


Fig.2: STEM ADF image (left) and atomic resolution image of the aluminum rich precipitate (right).



# Helium ion microscopy and sectioning of Spider Silk

Irina Iachina<sup>1,2</sup>, Jonathan R. Brewer<sup>2</sup>, Horst-Günter Rubahn<sup>1</sup>, Jacek Fiutowski<sup>1\*</sup>

<sup>1</sup>NanoSYD, Mads Clausen Institute, University of Southern Denmark

<sup>2</sup>Department of Biochemistry and Molecular Biology, University of Southern Denmark

\* corresponding author email: [fiutowski@mci.sdu.dk](mailto:fiutowski@mci.sdu.dk)

Focused ion beams have recently emerged as a powerful tool for ultrastructural imaging of biological samples. In this presentation, we will show that Helium Ion Microscopy (HIM), combined with ion milling, can be used to visualize the inner structure of both Major and Minor Ampullate silk fibres of the orb-web weaving spider *Nephila madagascariensis*. The internal nanofibrils were imaged in pristine silk fibres, with little or no observed damage to the sample structure. Furthermore, a method to cut/rupture the fibres using He<sup>+</sup> ions combined with internal sample tension is presented.

This work showed that HIM is a valuable tool for visualizing biological samples. The inherent sputtering of the He<sup>+</sup> ions could be used to mill the samples without damaging the biological structures softly. This contrasted with the Ne<sup>+</sup> ions, which caused damage to the spider silks' internal structure.

It was also shown that combining He<sup>+</sup> ion milling and the inherent tension in the spider silk sample made it possible to cut the specimen and visualize the rupturing process and the internal structures in the silk. This method could be used for other fibrous structures or even with non-elastic samples if used with stretched adhesive carbon tape, which could be used to apply the rupturing force.

The HIM images of the spider silk revealed that the rupturing process was highly dynamic involving rearrangement of the material in the fibre and showing strong indications of an internal fibril structure in the silk fibres with typical dimensions of 100-200nm. We anticipate that HIM will significantly contribute to some of the most challenging imaging applications and may open new directions in future bioimaging when paired with other imaging modalities.

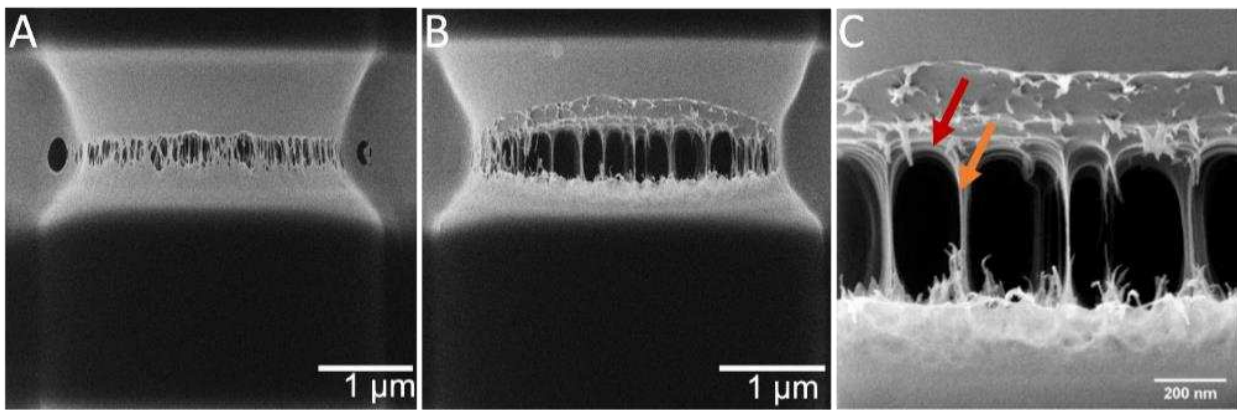


Fig. 1: HIM images of a MiS fiber that has been cut in half by He<sup>+</sup> ion sputtering. A) HIM image of a MiS fiber before breaking. B) image of the fiber stretching. C) Further zoom in on B) showing the individual fibrils being stretched (orange arrow), causing deformation in the base of the fibril (red arrow).

# Best Practices for Xe PFIB Preparation of Materials for Transmission Electron Microscopy

Suzy Vitale<sup>1</sup> and Joshua Sugar<sup>2</sup>

<sup>1</sup> Carnegie Institution for Science, Washington DC, 20015, USA

<sup>2</sup> Sandia National Laboratories, Livermore CA, 94550, USA

Over the past decade, there has been a concerted effort to develop gallium (Ga)-free focused ion beam (FIB) instruments for materials research. The helium ion microscope (HIM), neon ion microscope (NIM), and other similar inductively coupled plasma (ICP) multi-ion source focused ion beams such as oxygen (O) and argon (Ar) have been explored, though none have shown to be effective for the preparation of samples for transmission electron microscopy (TEM). Xenon (Xe) plasma focused ion beam (PFIB), thanks to its heavier ion mass (and thus high material sputter yield) and chemical inertness, has become the go-to instrument for large length scale (>100  $\mu\text{m}$ ) cross sectioning and FIB-SEM tomography. But can the Xe PFIB fabricate TEM samples appropriate for nanoscale microstructural characterization?

Xe PFIB has its own set of challenges relative to Ga FIB. Xe PFIB would not typically be the go-to tool for nanoscale milling applications like TEM sample preparation because the ICP Xe PFIB will produce a beam with a larger diameter than the liquid metal ion source (LMIS) Ga FIB in the range of currents typically used for TEM prep. A beam with a larger diameter can be more difficult to use in smaller scale (< 500 nm) applications. The Xe PFIB also has a greater range of current density across the width of the beam. The outer-most edges of the ion beam, called beam tails, are more prominent in the Xe PFIB compared to the Ga FIB whose beam tails are virtually non-existent in most modern FIB instruments. The Xe PFIB beam tails add complexity to finer scale milling applications such as thinning a TEM sample to electron transparency.

In this talk, we will demonstrate that it is possible to create Ga-quality TEM samples using Xe PFIB. This quality is defined as curtain-free and uniformly thin over any 25  $\mu\text{m}^2$  area where the relative thickness ( $t/\lambda$ ) measured with electron energy loss spectroscopy (EELS) is less than or equal to 1. We will identify best-practices and ideal parameters settings (including deposition energy (kV), deposition height, ion beam incident milling angle, and beam placement on the sample) to show that it is possible to fabricate two types of Xe PFIB-made samples: 1) a curtain-free uniformly thin electron transparent area where  $t/\lambda \approx 0.9$  over 38  $\mu\text{m}^2$ , and 2) a curtain-free non-uniformly thin electron transparent area where  $t/\lambda$  ranges from 0.3 to 1.1 over 48  $\mu\text{m}^2$ . Both sample types exhibit most, but not all, aspects of Ga-quality. Because of Xe PFIB's relatively large beam size and wide beam tails, a balance must be struck between ultimate sample thinness and the overall size of the electron transparent area.

## New Applications in Energy Research Enabled by a Triple Beam, Dual Chamber FIB with Isotropic Tomographic Voxels

Ben Tordoff, Ph.D.

EUFN 2022

Energy materials research relies heavily on a deep understanding of material and device microstructure to make advancements. Moreover, because these materials exist in complex devices constructed of many materials and interfaces between them, the critical information is often buried beneath the surface of a material specimen or encapsulated in a closed-form device. Furthermore, the final devices derive their performance from the nanoscale 3D arrangement and microstructures of the constituent materials involved. Recently, a commercial focused-ion beam scanning electron microscope was developed with an integrated fs-laser mill attached to the load lock of the microscope, opening the door to new and more powerful analysis approaches in energy materials research. Applications include rapid access to deeply buried structures for high resolution imaging and analytics, large area cross-section preparation, and massive material ablation for sample preparation of structures for, e.g. FIB-SEM tomography, TEM, APT, or X-ray nanoCT. Additionally, each of these may be correlatively guided by prior imaging with techniques like 3D X-ray microscopy to enable targeted analysis and preparation. By combining the laser mill with leading technologies for 3D FIB-SEM tomography with true isotropic voxels, this platform enables comprehensive quantitative analysis of energy materials at the nanoscale. We illustrate these concepts with a number of application examples and use cases demonstrating the utility of the approach.

# Application of FIB-TOF-SIMS for 3D high-resolution chemical characterization of Li-ion solid-state batteries

Agnieszka Priebe<sup>1\*</sup>, Abdessalem Aribia<sup>2</sup>, Jordi Sastre<sup>2</sup>, Moritz H. Futscher<sup>2</sup>, Yaroslav E. Romanyuk<sup>2</sup> and Johann Michler<sup>1</sup>

<sup>1</sup> Empa, Swiss Federal Laboratories for Materials Science and Technology, Laboratory for Mechanics of Materials and Nanostructures, Feuerwerkerstrasse 39, CH-3602 Thun, Switzerland

<sup>2</sup> Empa, Swiss Federal Laboratories for Materials Science and Technology, Laboratory for Thin Films and Photovoltaics, Überlandstrasse 129, Dübendorf, CH-8600 Switzerland

\* corresponding author email: [agnieszka.priebe@empa.ch](mailto:agnieszka.priebe@empa.ch)

The growing interest in Li-ion batteries (LIBs) is particularly driven by the vast need for electric vehicles (EVs), providing environment-friendly (i.e., renewable, clean and zero-emission) energy storage systems.[1] The development of new chemical systems, ensuring improved LIBs' capacity, power and efficiency as well as the reduction of production costs, was a great milestone. However, the safety issues as well as short lifetime are still the critical problems of the conventional LIBs. All-solid-state-batteries (ASSBs) using inorganic solid-state electrolytes (SSE) are great alternatives to commonly used LIBs based on highly flammable organic liquid electrolytes as they enable for operation in a wide range of temperatures.[2] Therefore, ASSBs are considered as one of the most promising future energy storage technologies. The characterization of ASSBs is demanding because they contain buried structures and heterogeneous interfaces. Furthermore, detection of Li and representing its 3D distribution with nanoscale resolution is attainable only to few analytical techniques. In this work, we exploit the outstanding potential of FIB-TOF-SIMS (focused ion beam time-of-flight secondary ion mass spectrometry, Figure 1) for comprehensive chemical characterization of novel Li-containing thin films, which are potential materials for the future generation batteries. This technique allows for parallel detection of all sample components with high spatial resolution (i.e., lateral resolution < 50 nm and depth resolution < 10 nm)[3–5] and high sensitivity (ppm). Our studies show that FIB-TOF-SIMS can reveal presence of 400±100 nm overlithiated grains and 100±30 nm nanoparticles with an increased  $^7\text{Li}^{16}\text{O}^+$  ion content in the Li- and Ni-rich layered oxide with the composition of  $\text{Li}_x\text{Ni}_{0.8}\text{Mn}_{0.1}\text{Co}_{0.1}\text{O}_2$  (LR-NMC811,  $x>1$ ) as well as monitor structural changes upon air exposure (both on the surface and in the bulk).[6] Furthermore, we demonstrate that simultaneous delivery of fluorine gas during FIB-TOF-SIMS can significantly improve the quality of acquired TOF-SIMS data.[7–11] The latter is proved using a novel Au/Li<sub>7</sub>La<sub>3</sub>Zr<sub>2</sub>O<sub>12</sub>/Pt/MgO/Si multilayer. In this case, the enhanced TOF-SIMS data helped understand the operation of the system. The detection of Li within Au layer after the polarization measurements explained the previously observed formation of internal electric field.[12] Our studies prove that FIB-TOF-SIMS is a powerful technique delivering essential insights into the complex structure of novel

Li-based materials, which can help optimize the functionality of future energy storage technologies.

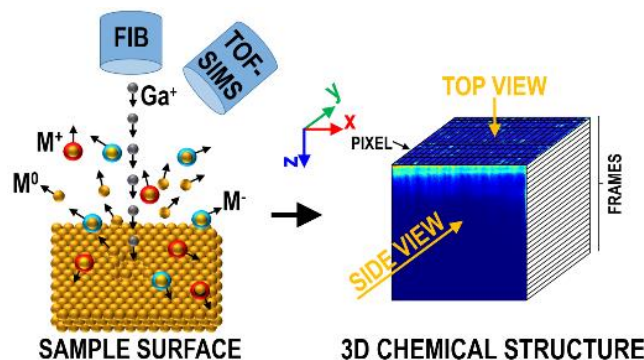


Fig. 1: FIB-TOF-SIMS is one of few analytical techniques allowing sample's 3D chemical structure to be assessed with nanoscale resolution.

- [1] Zeng X, Li M, Abd El-Hady D, Alshitari W, Al-Bogami A S, Lu J and Amine K 2019 Commercialization of Lithium Battery Technologies for Electric Vehicles *Adv. Energy Mater.* **9** 1–25
- [2] Zheng F, Kotobuki M, Song S, Lai M O and Lu L 2018 Review on solid electrolytes for all-solid-state lithium-ion batteries *J. Power Sources* **389** 198–213
- [3] Whitby J A, Östlund F, Horvath P, Gabureac M, Riesterer J L, Utke I, Hohl M, Sedláček L, Jiruše J, Friedli V, Bechelany M and Michler J 2012 High spatial resolution time-of-Flight secondary ion mass spectrometry for the masses: A novel orthogonal ToF FIB-SIMS instrument with in situ AFM *Adv. Mater. Sci. Eng.* **2012** 1–13
- [4] Priebe A, Barnes Jean-Paul, Edwards T E J, Huszár E, Pethö L and Michler J 2020 Elemental characterization of Al nanoparticles buried under a Cu thin film – TOF-SIMS vs. STEM/EDX *Anal. Chem.* **92** 12518–12527
- [5] Priebe A, Barnes J-P, Edwards T E J, Pethö L, Balogh I and Michler J 2019 3D imaging of nanoparticles in an inorganic matrix using TOF-SIMS validated with STEM and EDX *Anal. Chem.* **91** 11834–11839 Article
- [6] Priebe A, Aribia A, Sastre J, Romanyuk Y E and Michler J 2022 3D High-Resolution Chemical Characterization of Sputtered Li-rich NMC811 Thin Films Using TOF-SIMS *submitted*
- [7] Priebe A, Utke I, Pethö L and Michler J 2019 Application of a Gas-Injection System during the FIB-TOF-SIMS Analysis - Influence of Water Vapor and Fluorine Gas on Secondary Ion Signals and Sputtering Rates *Anal. Chem.* **91** 11712–22
- [8] Priebe A, Pethö L and Michler J 2020 Fluorine Gas Coinjection as a Solution for Enhancing Spatial Resolution of Time-of-Flight Secondary Ion Mass Spectrometry and Separating Mass Interference *Anal. Chem.* **92** 2121–9
- [9] Priebe A, Huszar E, Nowicki M, Pethö L and Michler J 2021 Mechanisms of fluorine-induced separation of mass interference during TOF-SIMS analysis *Anal. Chem.*
- [10] Priebe A, Xie T, Pethö L and Michler J 2020 Potential of gas-assisted time-of-flight secondary ion mass spectrometry for improving elemental characterization of complex metal-based systems *J. Anal. At. Spectrom.* **35** 2997–3006
- [11] Priebe A, Pethö L, Huszár E, Xie T, Utke I and Michler J 2021 High sensitivity of fluorine gas-assisted FIB-TOF-SIMS for chemical characterization of buried sublayers in thin films *ACS Appl. Mater. Interfaces*
- [12] Priebe A, Sastre J, Futscher M H, Jurczyk J, Puydinger Dos Santos M V., Romanyuk Y E and Michler J 2021 Detection of Au<sup>+</sup> Ions During Fluorine Gas-Assisted Time-of-Flight Secondary Ion Mass Spectrometry (TOF-SIMS) for the Complete Elemental Characterization of Microbatteries *ACS Appl. Mater. Interfaces*

# On demand spatially controlled fabrication of single photon emitters in Si

Gregor Hlawacek\*, Nico Klingner, Michael Hollenbach, Ulrich Kentsch, and Georgy Astakhov

<sup>1</sup> Institute for Ion Beam Physics and Materials Research, Helmholtz—Zentrum Dresden—Rossendorf, 01328 Dresden, Germany

\* corresponding author email: [g.hlawacek@hzdr.de](mailto:g.hlawacek@hzdr.de)

Single photon emitters (SPE) are fundamental building blocks for future quantum technology applications. However, many approaches lack the required spatial placement accuracy and Si technology compatibility required for many of the envisioned applications. Here, we present a method to fabricate at will placed single or few SPEs emitting in the telecom O-band in Silicon [1]. The successful integration of these telecom quantum emitters into photonic structures such as micro-resonators, nanopillars and photonic crystals with sub-micrometer precision paves the way toward a monolithic, all-silicon-based semiconductor-superconductor quantum circuit for which this work lays the foundations. To achieve our goal we employ home built AuSi liquid metal alloy ion sources (LMAIS) and an Orsay Physics CANION M31Z+ focused ion beam (FIB). Silicon-on-insulator substrates from different fabrication methods have been irradiated with a spot pattern. 6 to 500  $\text{Si}^{2+}$  ions have been implanted per spot using an energy of 40 keV. For the analysis and confirmation of the fabrication of true SPEs a home build photo luminescence setup has been used. G-centers formed by the combination of two carbon atoms and a silicon atom are confirmed by measurements of zero phonon lines (ZPL) at the expected wave length of 1278 nm for the case of carbon rich SOI wafers. In the case of ultra clean SOI wafers and high ion fluxes emission from tri-interstitial Si complexes is observed. The SPE nature of these so called W-centers has also been confirmed by ZPL measurements at 1218 nm. The achieved lateral SPE placement accuracy is below 100 nm in both cases and the success rate of SPE formation is more than 50%. After a discussion of the formation statistic we also present an approach how our FIB based approach can be upscaled to wafer-scale nanofabrication of telecom SPEs compatible with complementary metal oxide semiconductor (CMOS) technology for very large scale integration (VLSI).



[1] M. Hollenbach, N. Klingner, N. S. Jagtap, L. Bischoff, C. Fowley, U. Kentsch, G. Hlawacek, A. Erbe, N. V. Abrosimov, M. Helm, Y. Berencén, and G. V. Astakhov, "Wafer-scale nanofabrication of telecom single-photon emitters in silicon," (2022), arXiv:2204.13173 [quant-ph].

# Magnetic patterning using Ne, Co, and Dy FIB

K. Lenz<sup>1\*</sup>, J. Pablo-Navarro<sup>1</sup>, N. Kligner<sup>1</sup>, G. Hlawacek<sup>1</sup>, F. Samad<sup>1,2</sup>, R. Narkowicz<sup>1</sup>, R. Hübner<sup>1</sup>, A. Kákay<sup>1</sup>, H. Canzever<sup>1</sup>, W. Pilz<sup>3</sup>, F. Meyer<sup>3</sup>, P. Mazarov<sup>3</sup>, L. Bischoff<sup>1</sup>, R. Bali<sup>1</sup>, J. Lindner<sup>1</sup>

<sup>1</sup> Helmholtz-Zentrum Dresden-Rossendorf, Institute of Ion Beam Physics and Materials Research, Dresden, Germany

<sup>2</sup> Institute of Physics, University of Technology Chemnitz, Chemnitz, Germany

<sup>3</sup> Raith GmbH, Dortmund, Germany

\* corresponding author email: [k.lenz@hzdr.de](mailto:k.lenz@hzdr.de)

Magnetic nanostructures needed for magnonics and spintronics are usually processed by conventional lithography techniques in combination with lift-off or broad-beam ion etching. However, it has been shown [1] that the quality and shape of the structures' edges play an important role for the magnetization dynamics as structures become smaller and smaller. Furthermore, regarding optical measurement techniques, hard-to-remove resist masks that become hardened by ion etching are problematic. Direct-writing focused ion beams (FIB) do not have these issues. In addition, using non-standard ion species opens various paths for local magnetic patterning, i.e., influencing the magnetic properties locally.

I will present results for maskless magnetic patterning of ferromagnetic nanostructures using He and Ne ions as well as a few liquid metal alloy ion sources (LMAIS) for FIB systems. He/Ne FIBs are well established and commercially available. Irradiation of (paramagnetic) FeAl films by Ne ions creates local ferromagnetic nanostructures caused by disorder that are embedded in a paramagnetic matrix [2]. The precise Ne FIB also enables us to trim the edges of magnetic nanostructures enhancing their magnetic fidelity and creating certain localized magnon states at the edges of the samples. Using specifically developed LMAIS, like e.g., Co<sub>36</sub>Nd<sub>64</sub>, CoDy, or CuDy [3,4] in combination with a Wien mass filter offers further new paths for magnetic patterning. I will present results on the modification of Ni<sub>80</sub>Fe<sub>20</sub> (permalloy) strip samples. Using the CoNd LMAIS a narrow track of Co ions was implanted. The induced magnetic changes were measured with microresonator ferromagnetic resonance (FMR) before and after the implantation. Structures as small as 30 nm can be implanted up to a concentration of 10 % near the surface. Such lateral resolution is hard to reach for other lithographic methods. Using Dy ions one can locally increase the Gilbert damping parameter of the magnetization dynamics by more than a factor of four with a lateral resolution of about 100 nm.

[1] R. D. McMichael *et al.*, Phys. Rev. B 74, 024424 (2006)

[2] H. Cansever *et al.*, Sci. Rep., accepted (2022)

[3] L. Bischoff *et al.*, J. Vac. Sci. Technol. B, accepted (2022).

[4] L. Bischoff *et al.*, Nucl. Instrum. Method. Phys. Res. B 161-163, 1128 (2000).

# **fibTOF: The strength of SIMS capabilities on FIB-SEM microscopes**

Lex Pillatsch\*, James Whitby and Janos Metzger

TOFWERK AG, Schorenstrasse 39, CH-3645, Thun, Switzerland

\* corresponding author email: [lex.pillatsch@tofwerk.com](mailto:lex.pillatsch@tofwerk.com)

The combination of a FIB-SEM microscope with a mass analyzer provides a cost-effective route to high spatial resolution chemical imaging using secondary ion mass spectrometry (FIB-SIMS) [1]. The high mass resolving power, together with the ability to collect information about all elements simultaneously, has made the use of time-of-flight mass analyzers popular for this purpose, especially in the field of materials science. Because information about all mass to charge ratios is collected, retrospective analysis of the data sets can allow for new insights into local elemental distributions and correlations without repeating the measurement. FIB-SIMS is based on the direct measurement of the sputtered ions. The yields of secondary ions are in general high for positive ions from alkali metals and, also for negative ions from halogens. This is especially true for lithium and fluorine, which can be difficult to map by other techniques. FIB-SIMS is therefore a technique well placed to support research of current and next-generation rechargeable battery materials (e.g., for lithium ion, sodium ion and fluoride ion batteries).


The correct identification of the sputtered ions is mandatory for the characterization of the material properties. Knowledge about the possible sample composition is required for the data interpretation. A high mass resolving power is essential to resolve molecules from elements at the same atomic mass unit. Molecular fragments containing the element of interest can confirm the correct elemental identification.

FIB-SIMS is thus a powerful technique that can complement other characterization methods (electrochemical, XRD, SEM imaging, or XPS). Having the specimen within the vacuum chamber of a FIB-SEM microscope means that a fresh (not oxidized) surface can be prepared, although the use of a cryostage and grazing incidence FIB beam may be necessary to get the best quality surface when using Ga<sup>+</sup> primary (FIB) ions [3].

Recent results from rechargeable battery materials recorded by the fibTOF and methods for the correct data interpretation will be presented.

## References

- [1] L. Pillatsch, F. Östlund, J. Michler "FIBSIMS: A review of secondary ion mass spectrometry for analytical dual beam focused ion beam instruments" *Progress in Crystal Growth and Characterization of Materials*, Volume 65, Issue 1, February 2019, Pages 1-19. <https://doi.org/10.1016/j.pcrysgrow.2018.10.001>
- [2] J. Sastre, M. H. Futscher, L. Pompizi, A. Aribia, A. Priebe, J. Overbeck, M. Stiefel, A. N. Tiwari, Y. E. Romanyuk "Blocking lithium dendrite growth in solid-state batteries with an ultrathin amorphous Li-La-Zr-O solid electrolyte" *Communications Materials*, Vol. 2, pp. 1–10 (2021) <https://www.nature.com/articles/s43246-021-00177-4>
- [3] Lee, J. Z. et al. "Cryogenic focused ion beam characterization of lithium metal anodes". *ACS Energy Lett.* 4, 489–493 (2019). <https://doi.org/10.1021/acsenergylett.8b02381>
- [4] Ch. Jiao, L. Pillatsch, J. Mulders, D. Wall "Three-Dimensional Time-of-Flight Secondary Ion Mass Spectrometry and DualBeam FIB/SEM Imaging of Lithium-ion Battery Cathode", *Microscopy and Microanalysis*, Vol. 25(S2), pp. 876-877 (2019) <https://doi.org/10.1017/S1431927619005117>

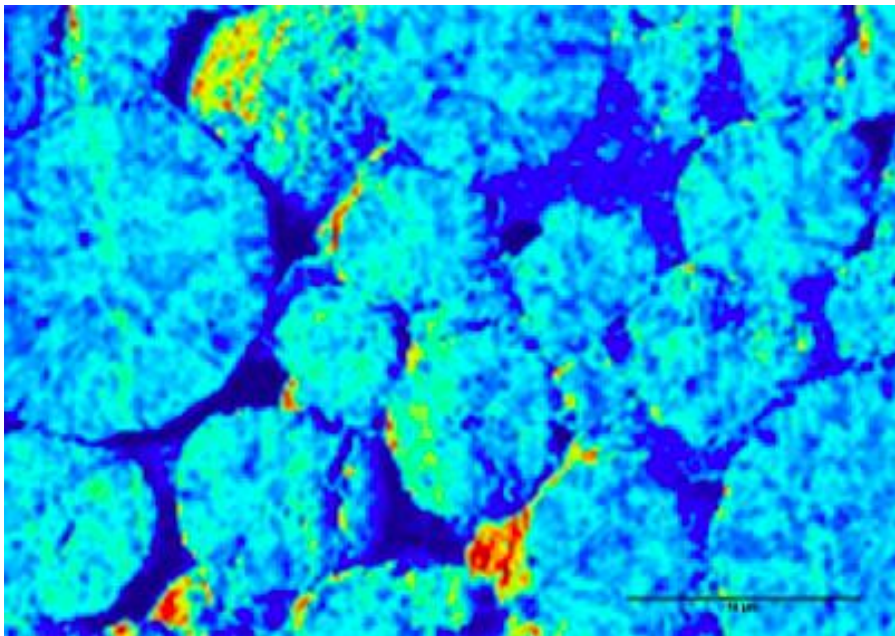

*Figure 1: Lithium distribution in a  $\text{LiNi}_x\text{Mn}_y\text{Co}_{1-x-y}\text{O}_2$  (NMC) cathode, from [4].*

# Positioned generation of luminescence defects in 2D materials by helium ion beams

A. Hötger,<sup>1</sup> T. Amit,<sup>2</sup> J. Klein,<sup>3</sup> K. Barthelmi,<sup>1</sup> T. Pelini,<sup>4</sup> A. Delhomme,<sup>4</sup> S. Rey,<sup>5</sup> M. Potemski,<sup>4,6</sup> C. Faugeras,<sup>4</sup> G. Cohen,<sup>2</sup> D. Hernangómez-Pérez,<sup>2</sup> T. Taniguchi,<sup>7</sup> K. Watanabe,<sup>8</sup> C. Kastl,<sup>1</sup> S. Refaely-Abramson,<sup>2</sup> A.V. Stier,<sup>1</sup> J.J. Finley,<sup>1</sup> and A.W. Holleitner<sup>1\*</sup>

<sup>1</sup> Walter Schottky Institute and Physics Department, TU Munich, 85748 Garching, Germany.

<sup>2</sup> Department of Molecular Chemistry and Materials Science, Weizmann Institute of Science, Rehovot, Israel.

<sup>3</sup> Department of Materials Science and Engineering, Massachusetts Institute of Technology, Cambridge, Massachusetts 02139, USA

<sup>4</sup> Laboratoire National des Champs Magnétiques Intenses, CNRS-UGA-UPS-INS-EMFL, 38042 Grenoble, France.

<sup>5</sup> Department of Photonics Engineering, Technical University of Denmark, 2800 Kgs. Lyngby, Denmark.

<sup>6</sup> Institute of Experimental Physics, Faculty of Physics, University of Warsaw, 02-093 Warszawa, Poland

<sup>7</sup> International Center for Materials Nanoarchitectonics, National Institute for Materials Science, 1-1 Namiki, Tsukuba 305-0044, Japan.

<sup>8</sup> Research Center for Functional Materials, National Institute for Materials Science, 1-1 Namiki, Tsukuba 305-0044, Japan.

\* corresponding author email: [holleitner@wsi.tum.de](mailto:holleitner@wsi.tum.de)

Single spin defects in 2D transition-metal dichalcogenides are natural spin-photon interfaces for quantum applications. The creation of such point-defects in bulk crystals lacks either lateral or vertical precision. We overcome this disadvantage by irradiating atomically thin materials by a focused He-ion beam. Here we report on deterministic generation of optically active defects with a helium ion microscope in monolayer MoS<sub>2</sub>. In photoluminescence (PL) measurements on the irradiated sites we measure sharp emission lines ~200 meV below the optical bandgap of MoS<sub>2</sub>. [1] The He-ion beam dose can be decreased to create single emission lines at each irradiated spot with a yield of ~20%. These single emission lines emit single photons, which could be unambiguously proven in second-order correlation measurements. [2] In high-field magneto-photoluminescence we attribute the emission lines to single sulfur vacancies by combining experiment and ab-initio calculations. In the experiment, we reveal the lifting of spin-degeneracy of the involved defect bands even at zero magnetic field. These results highlight that defects in 2D semiconductors may be utilized for quantum technologies. [3]

[1] J. Klein et al.; Site-selectively generated photon emitters in monolayer MoS<sub>2</sub> via local helium ion irradiation, *Nature Comm.* 10 (2019), 2755

[2] J. Klein, L. Sigl et al.; Engineering the Luminescence and Generation of Individual Defect Emitters in Atomically Thin MoS<sub>2</sub>; *ACS Photonics* 8 (2021), 2

[3] A. Hötger et al.; Spin-defect characteristics of single sulfur vacancies in monolayer MoS<sub>2</sub>; *arXiv:2205.10286* (2022)

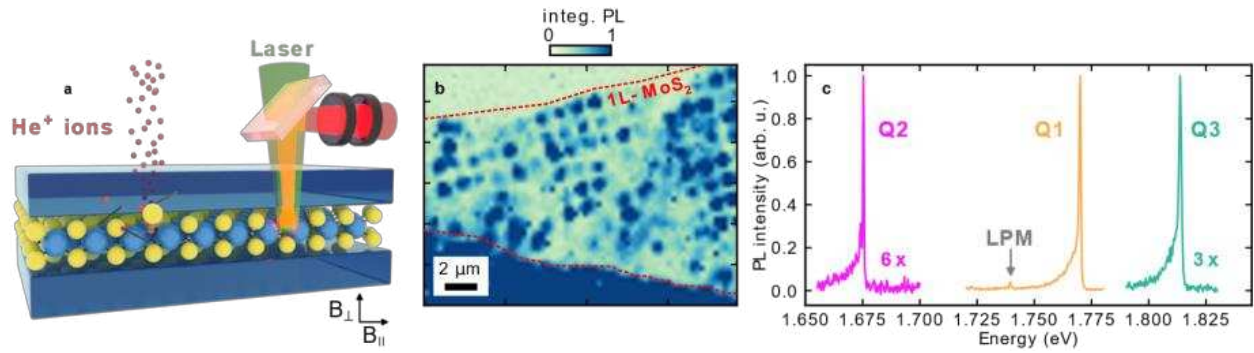


Fig. 1: a) Sketch of the He-ion irradiation and the magneto-optical measurement on the generated defect center. b) False color map of the defect luminescence showing the irradiated array pattern. c) Typical photoluminescence spectra of three different emission bands showing similar asymmetric line shapes.

# Fabrication of microstructured devices for grain boundary investigations in unconventional superconductor CeCoIn<sub>5</sub>

S. Mishra<sup>\*</sup>, S. M. Thomas, R. McCabe, S-Z Lin, E. D. Bauer and F. Ronning

<sup>1</sup> Los Alamos National Laboratory, Los Alamos, New Mexico, 87545, USA

<sup>\*</sup> corresponding author email: [sanu@lanl.gov](mailto:sanu@lanl.gov)

The study of grain boundaries (GB) in superconductors has both fundamental and applied interests. In high-temperature cuprate superconductors studies of the critical currents ( $J_c$ ) across GBs have provided important information on the symmetry of the superconducting order parameter and are critical for the observation of spontaneously generated half-flux magnetic quanta [1,2]. Similar to cuprate superconductors, heavy fermion superconductors (HFS) host rich physics in the form of unconventional superconducting phases with nodal quasiparticles. However, there have been relatively few phase-sensitive measurements of the superconducting order parameter thereby emphasizing the need for investigations of  $J_c$  across GBs in HFS.

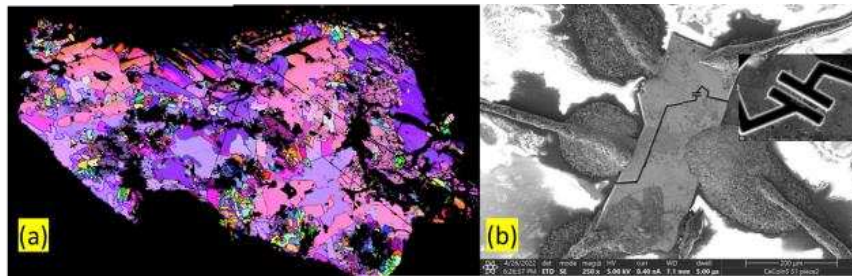


Fig. 1: (a) Shows an EBSD image of the differently oriented grains present in polycrystalline CeCoIn<sub>5</sub> sample and (b) shows the SEM image of a microstructured fabricated across a 90° GB using FIB milling.

In this talk, I will present results on GBs in polycrystalline samples of the HFS CeCoIn<sub>5</sub>. Electron backscatter diffraction (EBSD) performed on well-polished samples of polycrystalline CeCoIn<sub>5</sub> reveal that majority of grains are not randomly oriented as one would expect but grow at a misorientation angle of 90° with respect to their neighboring grain. We performed  $J_c$  studies across various such GBs by fabricating microstructured devices using focused ion-beam milling. Our investigations are crucial in understanding the superconducting order parameter symmetry of CeCoIn<sub>5</sub> and its potential use in devices for quantum information science.

[1] H. Hilgenkamp et al., Rev. Mod. Phys. 74, 485, 2002.

[2] C. C. Tsuei and J. R. Kirtley, Rev. Mod. Phys. 72, 969, 2000.





## ***Exploring Layered Conductors by 3D FIB micro-machining.***


*Carsten Putzke*

Some of the most enigmatic correlated states in the field of quantum matter arise in lower dimensions such as quasi-2d and quasi-1d materials. The layered crystal structure in organic and high- $T_c$  superconductors is an example of these materials, where a quasi-2D electronic structure gives rise to several unconventional electronic instabilities. The layered crystal structure allows for exfoliation in some materials which helped to uncover electrical transport properties such as the quantum Hall effect as well as enables spectroscopic probes in single layer and few layer systems. It is the same crystallographic anisotropy that hinders the study of interlayer electrical transport and spectroscopic probes such as angle resolved photoemission spectroscopy (ARPES).

In my talk I will demonstrate the novel experimental capabilities that focused ion beam micro-structuring enables in studying layered conductors. Confining the in-plane dimension of quasi-2D high purity metals of the Delafossites to length scales smaller than the electron mean free path gives rise to a novel realization of the particle-wave duality<sup>1,2</sup>.

Beyond these exciting physical phenomena of finite size confined pillars, novel experimental possibilities will be presented which enable previously inaccessible insight into layered materials.

1. M.D. Bachmann, et al. Nature Physics (accepted 2022), arXiv:2103.01332
2. C. Putzke, et al. Science 368, 6496 (2020)

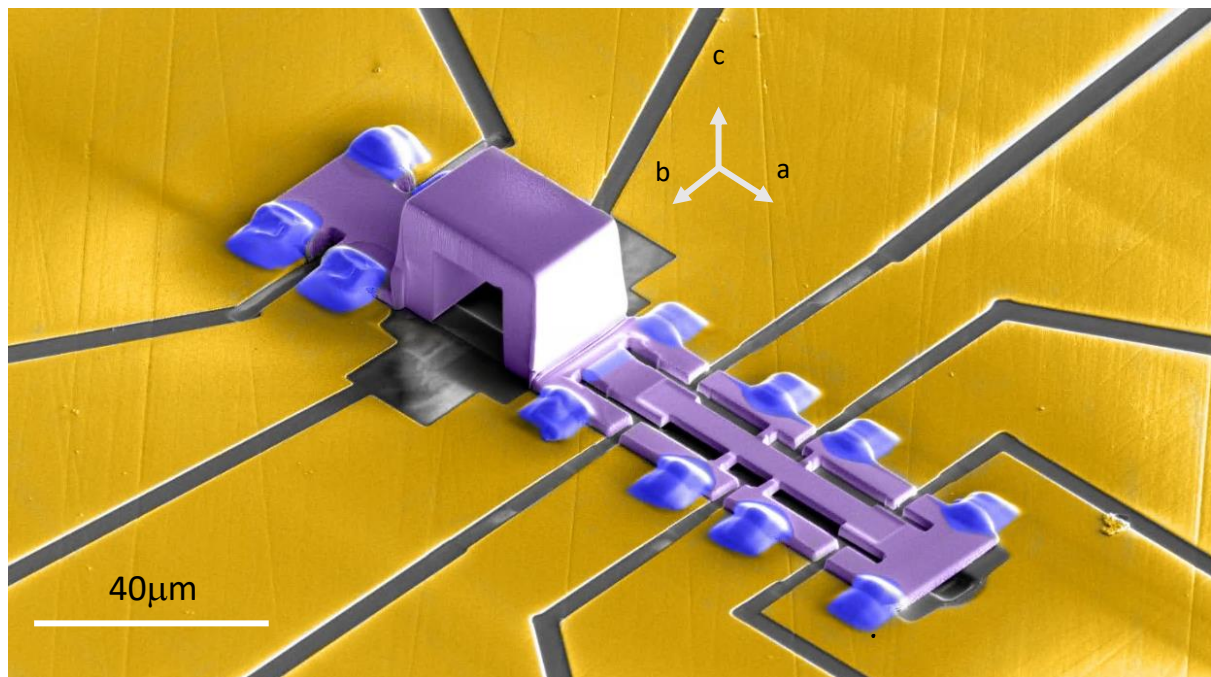

Scanning electron microscope image of a 3D PdCoO<sub>2</sub> micro-structure enabling to study the interplay between inter- and intra-layer electrical transport.

# Investigation of the Interaction of a Ga<sup>+</sup> Focused Ion Beam with Zirconia by Electron Backscatter Diffraction

N. Brachhold<sup>1\*</sup>, C.G. Aneziris<sup>1</sup>

<sup>1</sup> TU Bergakademie Freiberg, Institute of Ceramics, Refractories and Composite Materials, 09599 Freiberg, Germany

\* corresponding author email: [nora.brachhold@ikfww.tu-freiberg.de](mailto:nora.brachhold@ikfww.tu-freiberg.de)

The interaction of a Ga<sup>+</sup> focused ion beam with zirconia materials (Mg-PSZ, Y-PSZ) was investigated with the aim to identify suitable parameters for the sample preparation for surface sensitive analytical methods such as EBSD. The FIB parameters glancing angle (angle between ion beam and surface) and applied power (defined by acceleration voltage and ion current) were varied between two levels.

It was observed that the phase stability of monoclinic zirconia depended strongly on the FIB treatment parameters [1-2]. For a glancing angle of 5° and an applied acceleration voltage of 30 kV with an ion current of 30 nA a transformation of originally monoclinic grains to the tetragonal/cubic phase took place (Fig. 1). This was related to the implantation of Ga<sup>+</sup> ions and their presence as Ga<sub>2</sub>O<sub>3</sub> which works as stabilizer for the high temperature phases of zirconia [3-4].

Furthermore, the statistical analysis yielded that the influence of the FIB parameters on the surface quality corresponded well to the literature. The data showed that mild parameter levels (lower glancing angle and lower power of the ion beam, such as 5° and acceleration voltage of 5 kV combined with ion current of 4.8 nA, respectively) reduced the surface damage.

[1] H. Berek, C.G. Aneziris; *Effect of focused ion beam sample preparation on the phase composition of zirconia*; Ceramics International 44 (2018), 176435.

[2] N. Brachhold, H. Berek, J. Fruhstorfer, C.G. Aneziris; *Focused Ion Beam Parameters for the Preparation of Oxidic Ceramic Materials*; Advanced Engineering Materials 23 (2021), 2001235.

[3] A. Surpi, E. Göthelid, T. Kubart, D. Martin, J. Jensen; *Localised modifications of anatase TiO<sub>2</sub> thin films by a Focused Ion Beam*; Nuclear Instruments and Methods in Physics Research, Section B: Beam Interactions with Materials and Atoms 268 (2010), 3142.

[4] T. Yamamoto, A. Kurimoto; *Ga Ion-doped ZrO<sub>2</sub> Catalyst Characterized By XRD, XAFS, and 2-butanol decomposition*; Analytical Sciences 36 (2020), 41.

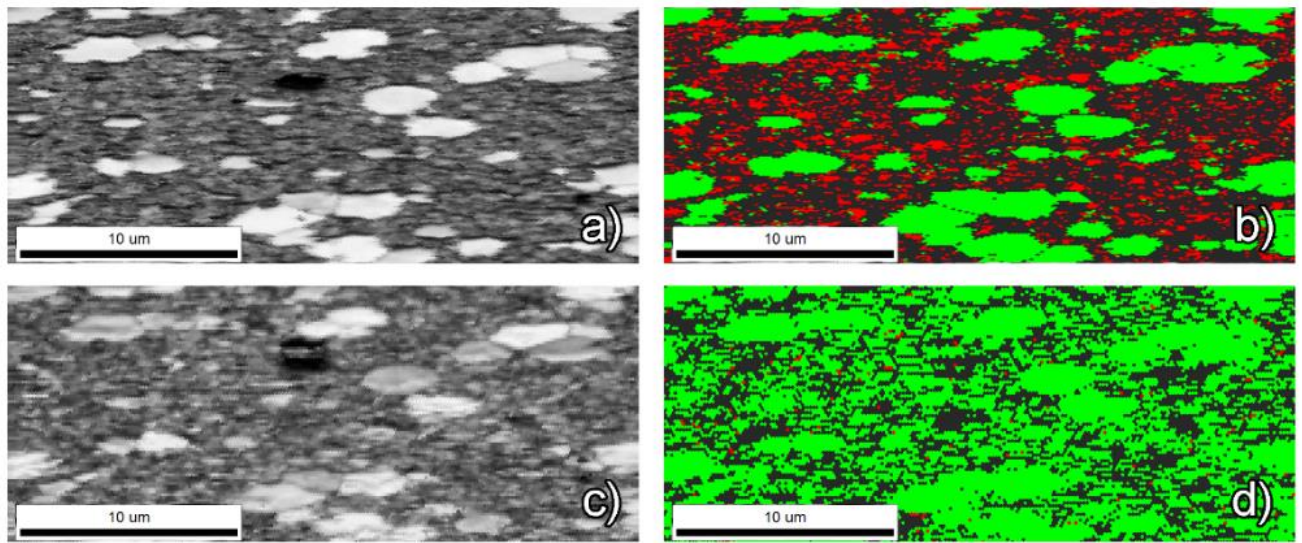


Fig. 1: Image quality (a, c) and phase distribution (b, d – red: monoclinic  $\text{ZrO}_2$ , green: cubic/tetragonal  $\text{ZrO}_2$ ) of Y-PSZ, (a,b) initial state, (c,d) state after a  $\text{Ga}^+$  FIB treatment at 30 kV and 30 nA with a glancing angle of  $5^\circ$